


# AUXILOCK® ADVANCE SYNDESMOSIS ULTIMATE KIT

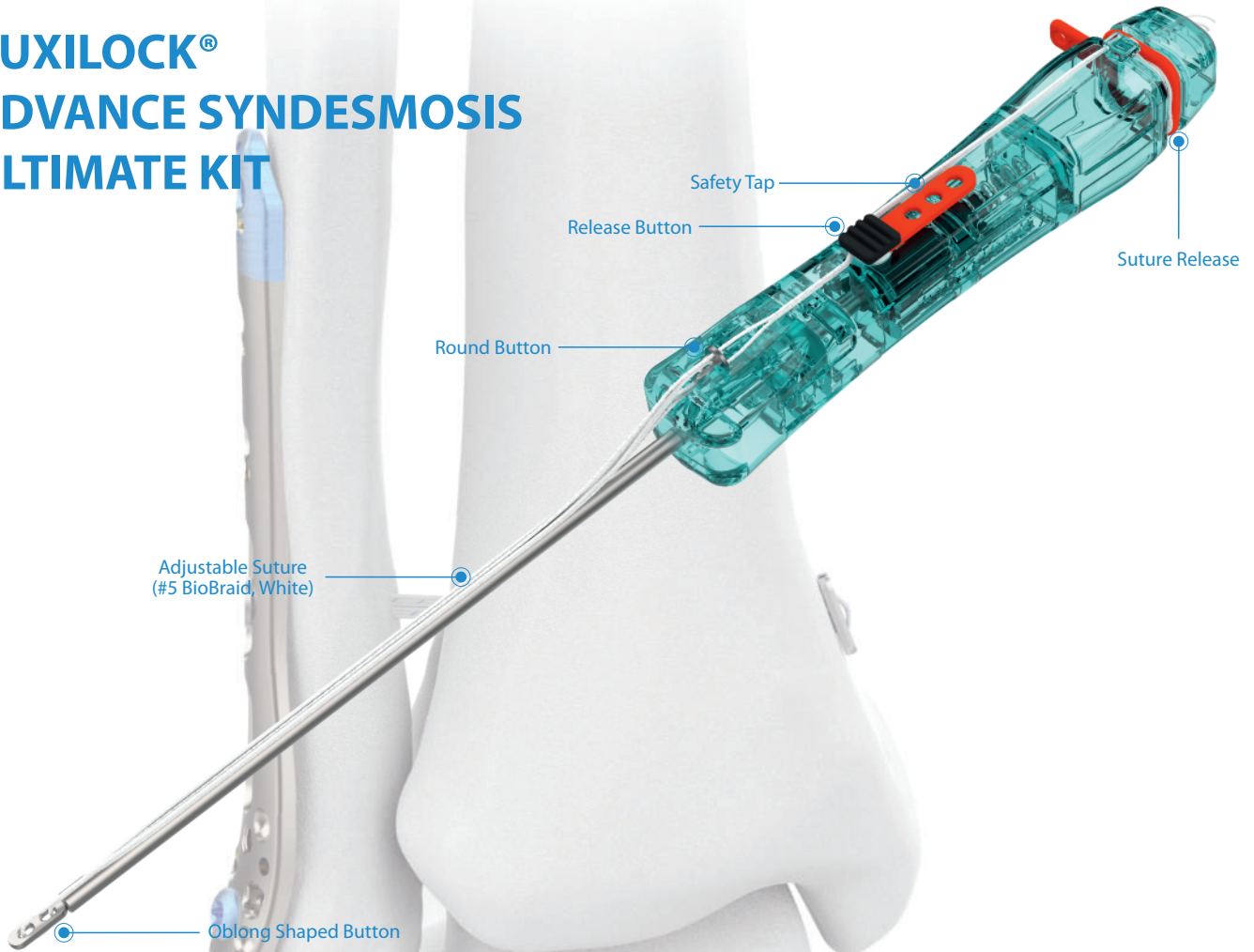

The Advance Syndesmosis Handle implant system features a unique delivery mechanism that allows the user to insert a Syndesmosis implant without pulling a needle through the medial skin. Advance Syndesmosis Handle system gives the user control to flip the medial oblong button below the skin for less soft-tissue interposition.

Ergonomic Advance Syndesmosis Handle implant tensioning handles have also been added to the system. These handles are designed to tension the Advance Syndesmosis Handle implant without hurting the surgeon's fingers. The Advance Syndesmosis Handle implant system also features a new trocar style drill bit that reduces skiving off the fibula when drilling in the transmalleolar plane 30° anterior to the coronal plane.

## Ordering Information

Code	Product Description
6-051-01	AUXILOCK® Advance Syndesmosis Ultimate Kit

## Related Instrumentation (Include in Kit)

Code	Product Description
7-147-01	Cannulated Drill Bit with Quick Coupling End, Ø3.5mm x Length 175mm
7-147-02	Ø1.5mm K-Wire

**STERILE** EO