


AUXILOCK® GFS ULTIMATE LARGE

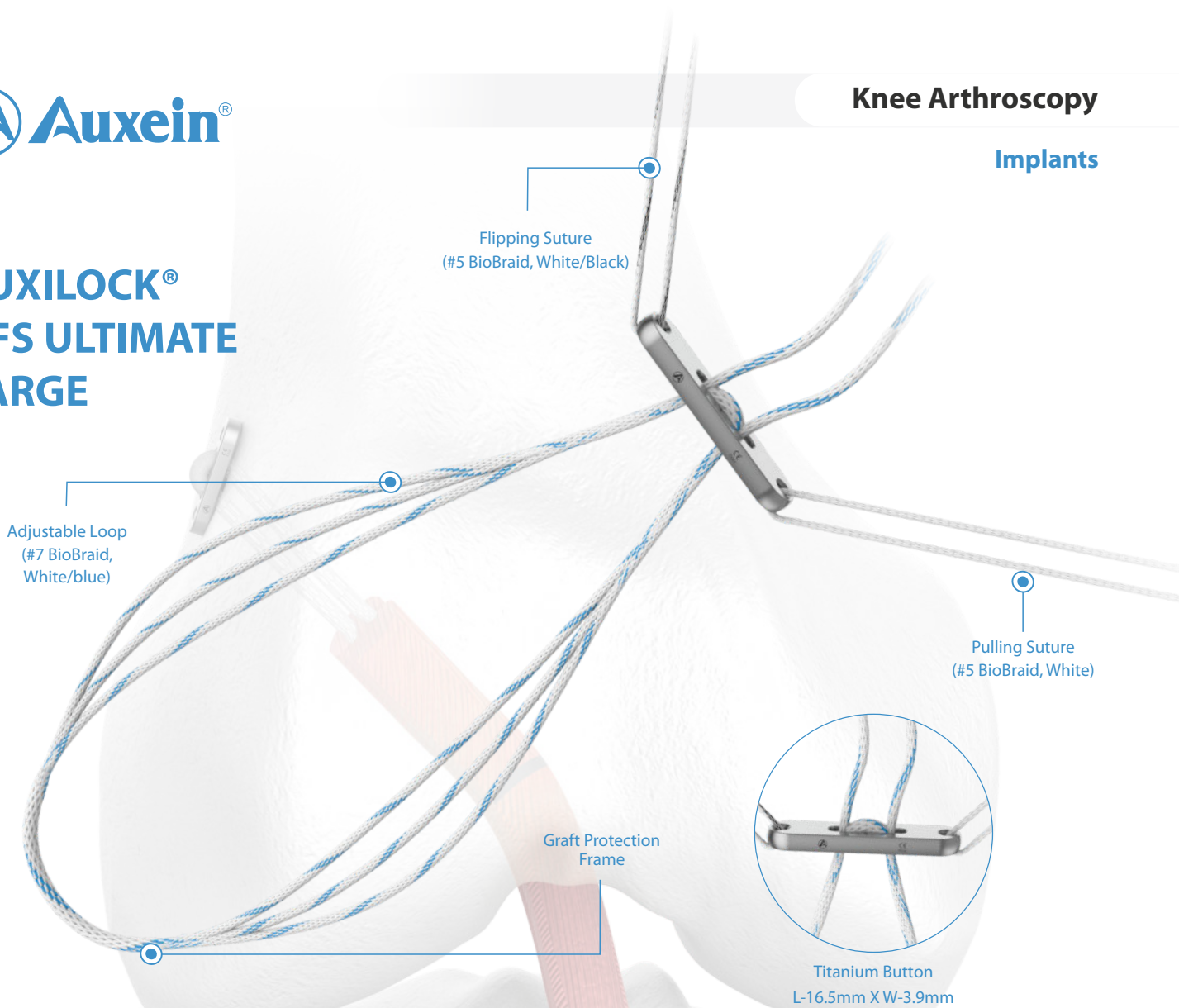

GFS Ultimate Large is an adjustable loop with an oblong shaped titanium button used for Cruciate Graft Reconstructions. GFS Ultimate provides a double locking mechanism which eliminates the need for knot tying. Hence, GFS Ultimate Large is considered as an excellent choice for revision surgeries.

GFS Ultimate Large provides three sutures (UHMWPE):

- Adjustable Suture (#7, White/Blue): It allows the surgeon to maximize the amount of graft inside the femoral tunnel, thereby optimizing the healing process. It also enables calibration of the loop to its optimum size.
- Pulling Suture (#5, White): It is available to pull the graft inside the tunnel.
- Flipping Suture (#5, White/Black): It ensures the flipping of the button.

Ordering Information

Code	Product Description	Length of Loop (mm)
6-006-02	AUXILOCK® GFS Ultimate Large Button, Button L: 16.5mm, W: 3.9mm (Adjustable Loop)	Max: 70, Min: 12

Recommended Instrumentation

Code	Product Description
7-013-57	Ø4.5mm Cannulated Headed Reamer

STERILE EO