


AUXILOCK® KNOTLESS PEEK OPTIMA PUSH-IN SUTURE ANCHOR

Shoulder Instability Repair

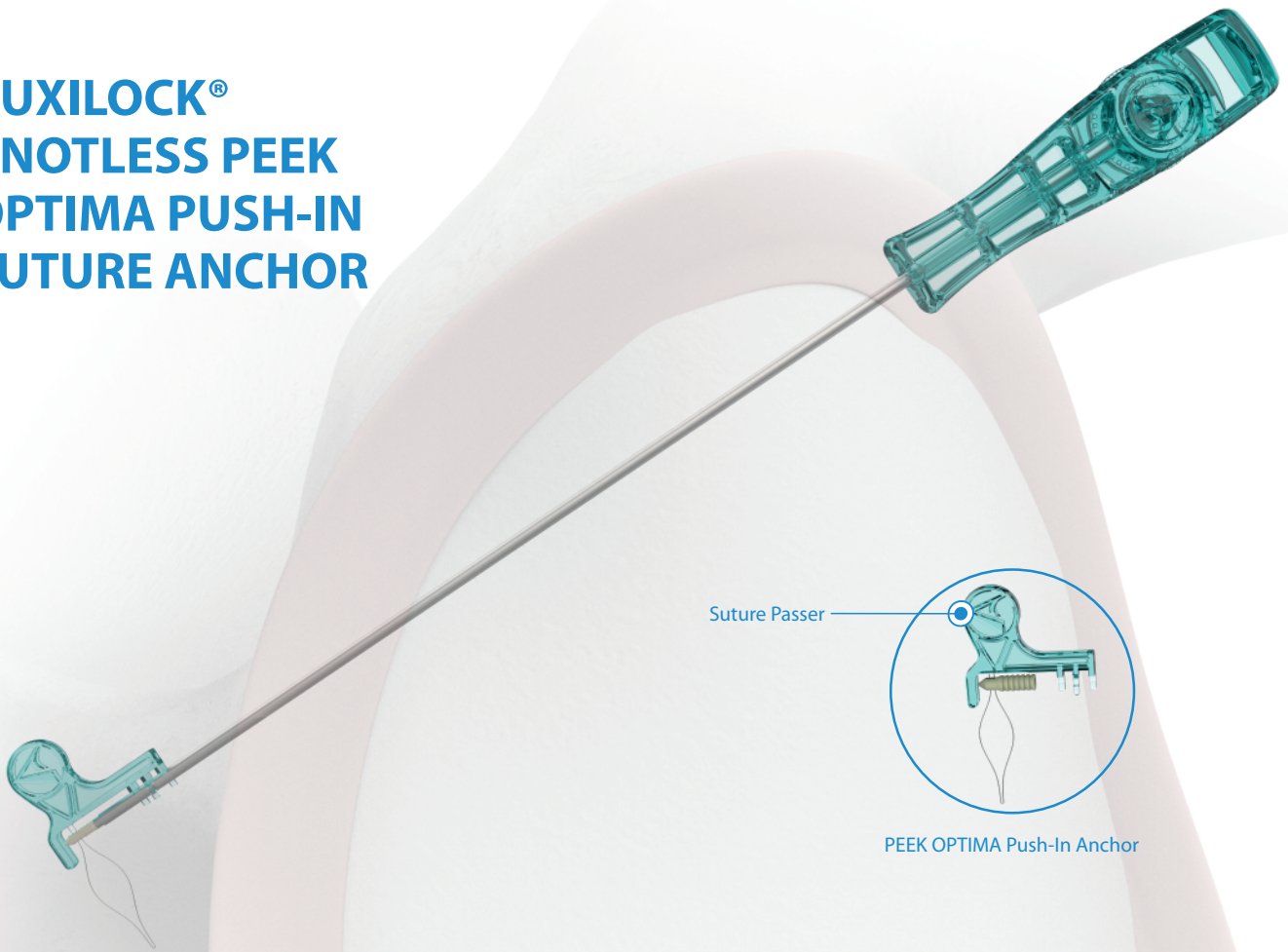

AUXILOCK® knotless push-in anchors are a step-saving alternative to conventional knotted suture anchors and are recommended for use in Bankart and SLAP repairs. The PEEK OPTIMA Push-In Suture Anchors are made of PEEK-OPTIMA (Poly-ether-ether-ketone). These Anchors are available in a variety of diameters of 2.8mm and 3.5mm providing intraoperative flexibility. The PEEK OPTIMA Push-In Suture Anchor is a one-piece body anchor with a material eyelet. It provides superior abrasion resistance due to PEEK's low coefficient friction. It is Biocompatible, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation.

Ordering Information

Code	Product Description
6-031-01	AUXILOCK® 2.8mm Knotless PEEK OPTIMA Push-In Suture Anchor With Driver, With Suture Passer
6-031-02	AUXILOCK® 3.5mm Knotless PEEK OPTIMA Push-In Suture Anchor With Driver, With Suture Passer

Recommended Instrumentation

Code	Product Description	Compatible Implants
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-031-01, 6-031-02
7-013-54	Ø2.5mm X 200 mm Drill Bit (For Knotless Anchor)	6-031-01
7-013-55	Ø3.2mm X 200 mm Drill Bit (For Knotless Anchor)	6-031-02
7-013-02	8 Point Drill Guide, Length 100mm	6-031-01, 6-031-02

STERILE EO