




AUXILOCK® ADVANCE SYNDESMOSIS ULTIMATE KIT

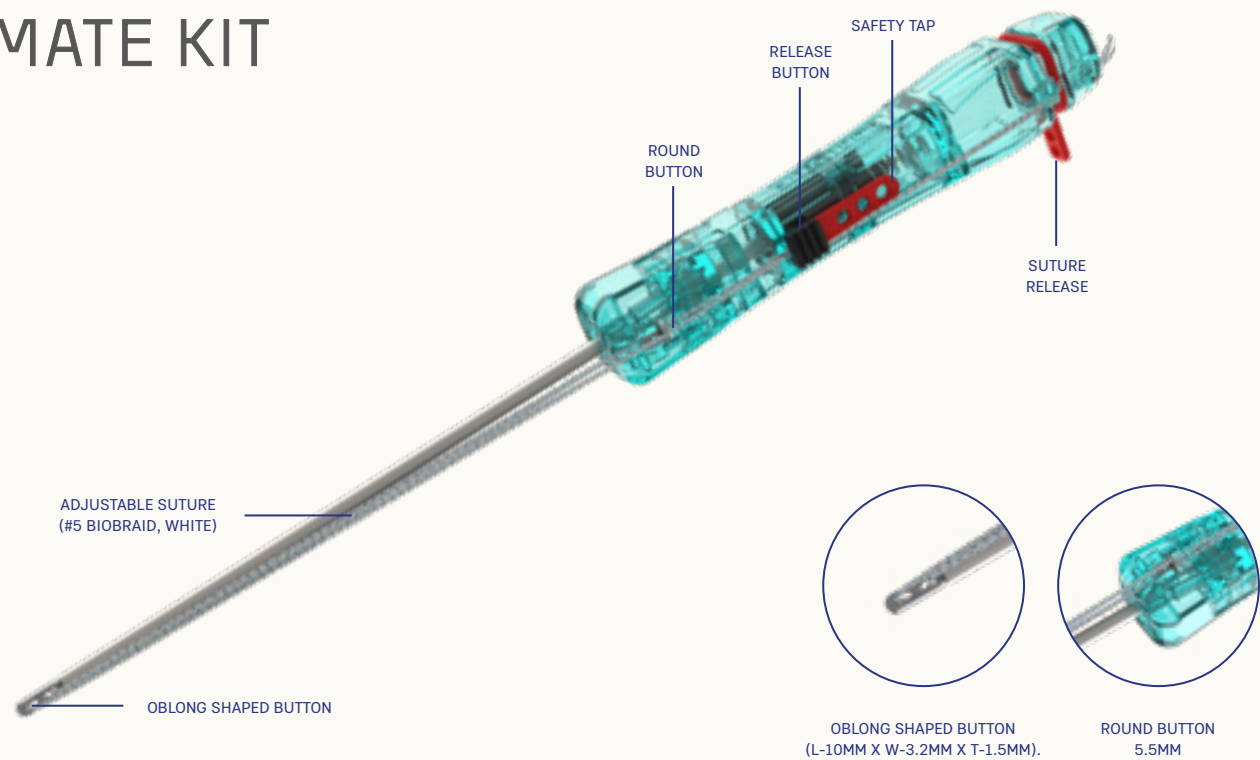

The AUXILOCK® Advance Syndesmosis implant system features a unique delivery mechanism that allows the user to insert a Syndesmosis implant without pulling a needle through the medial skin. Advance Syndesmosis Handle system gives the user control to flip the medial oblong button below the skin for less soft-tissue interposition.

Ergonomic Advance Syndesmosis implant tensioning handles have also been added to the system. The handle is designed tension the syndesmosis implant without hurting the surgeon fingers.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-051-01	AUXILOCK® Advance Syndesmosis Ultimate Kit

RELATED INSTRUMENTATION (INCLUDE IN KIT)

CODE	PRODUCT DESCRIPTION
7-147-01	Cannulated Drill Bit with Quick Coupling End, Ø3.5mm x Length 175mm
7-147-02	Ø1.5mm K-Wire

FOOT AND ANKLE