




AUXILOCK® AUXFIX BUTTON

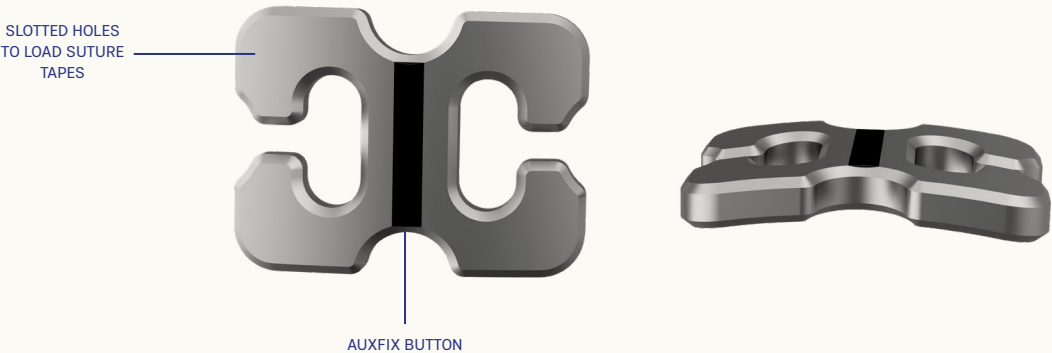

AUXILOCK® AuxFix Button is a titanium button with holes that allow use of multiple suture tapes for AC joint reduction. It provides very strong compression on the clavicles. The button is laid on the outer surface of the clavicle and underneath the coracoid. Only the suture material is passed through the clavicle and coracoid tunnel ensuring minimal bone loss during the procedures. Designed for use in chronic AC joint cases.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-037-01	AUXILOCK® AuxFix Button (AC Joint)

RELATED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-014-68	Passing Pin Ø2.4mm X 250mm
7-013-57	Ø4.5mm Cannulated Headed Reamer

AC JOINT REPAIR