




AUXILOCK® GFS MINI

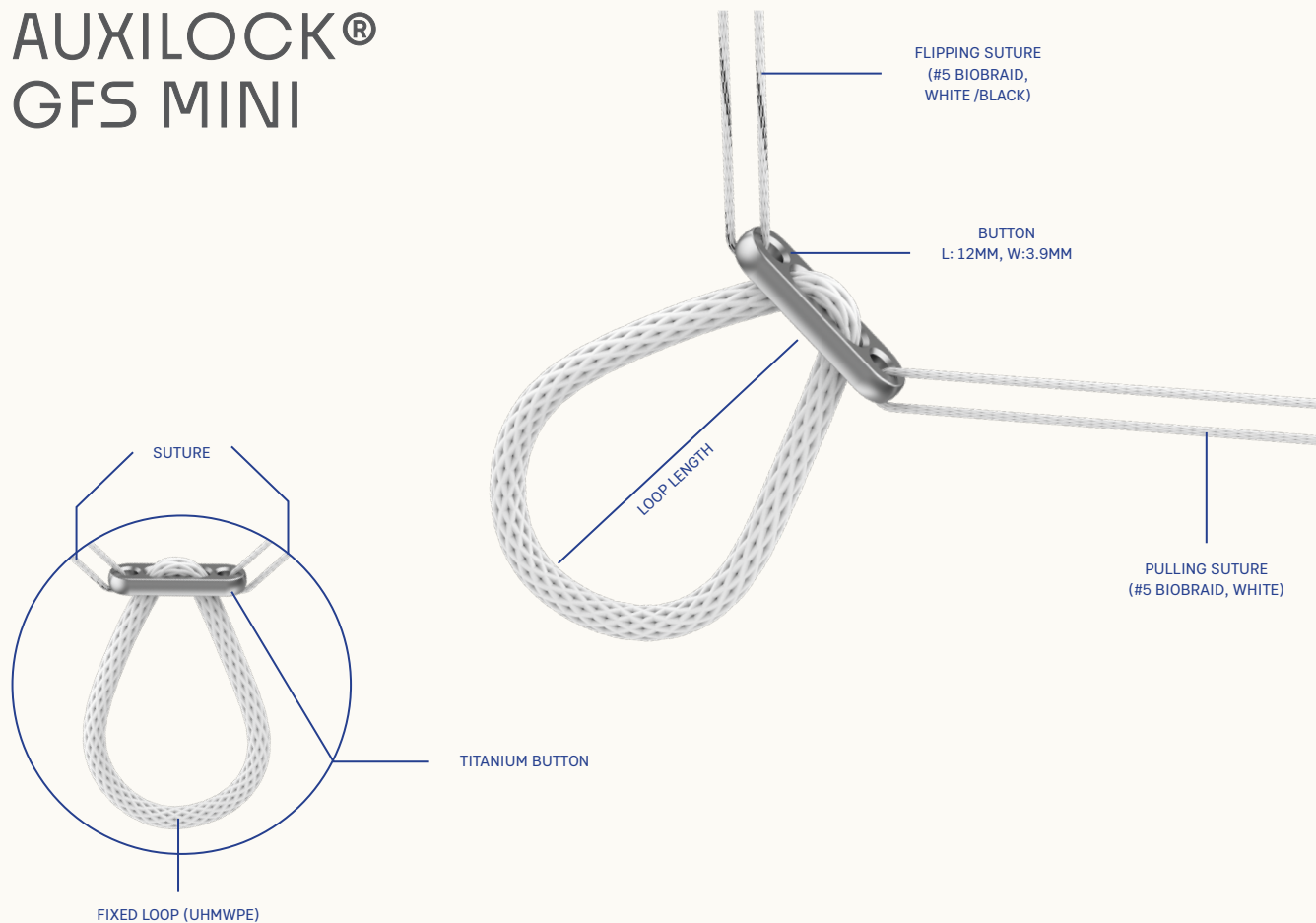

GFS Mini is a fixed loop with an oblong shaped button used for Graft Fixations in ACL/PCL Reconstruction. Being a continuous loop without any joint, GFS Mini provides a stronger fixation while eliminating the need for knot tying

The pre-loaded pulling and flipping braided sutures are available to ensure controlled pulling and flipping of the button in the transosseous tunnel. A 4.5mm cannulated headed reamer is provided in the instrument set for drilling the tunnel, allowing the GFS Mini button to pass easily into it. The fixed loop is made of UHMWPE (Ultra High Molecular Weight Polyethylene) and the oblong Button is made of titanium material.

- **Pulling Suture (#5, White):** It is available to pull the graft inside the tunnel.
- **Flipping Suture (#5, White/Black):** It ensures the flipping of the button.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-007-01	AUXILOCK® GFS Mini, Loop: 12mm, Button L: 12mm, W: 3.9mm
6-007-02	AUXILOCK® GFS Mini, Loop: 15mm, Button L: 12mm, W: 3.9mm
6-007-03	AUXILOCK® GFS Mini, Loop: 20mm, Button L: 12mm, W: 3.9mm
6-007-04	AUXILOCK® GFS Mini, Loop: 25mm, Button L: 12mm, W: 3.9mm
6-007-05	AUXILOCK® GFS Mini, Loop: 30mm, Button L: 12mm, W: 3.9mm
6-007-06	AUXILOCK® GFS Mini, Loop: 35mm, Button L: 12mm, W: 3.9mm
6-007-07	AUXILOCK® GFS Mini, Loop: 40mm, Button L: 12mm, W: 3.9mm
6-007-08	AUXILOCK® GFS Mini, Loop: 45mm, Button L: 12mm, W: 3.9mm
6-007-09	AUXILOCK® GFS Mini, Loop: 50mm, Button L: 12mm, W: 3.9mm
6-007-10	AUXILOCK® GFS Mini, Loop: 55mm, Button L: 12mm, W: 3.9mm

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-57	4.5mm Cannulated Headed Reamer	6-007-01, 6-007-02, 6-007-03, 6-007-04, 6-007-05, 6-007-06, 6-007-07, 6-007-08, 6-007-09, 6-007-10