




AUXILOCK® GFS ULTIMATE LARGE

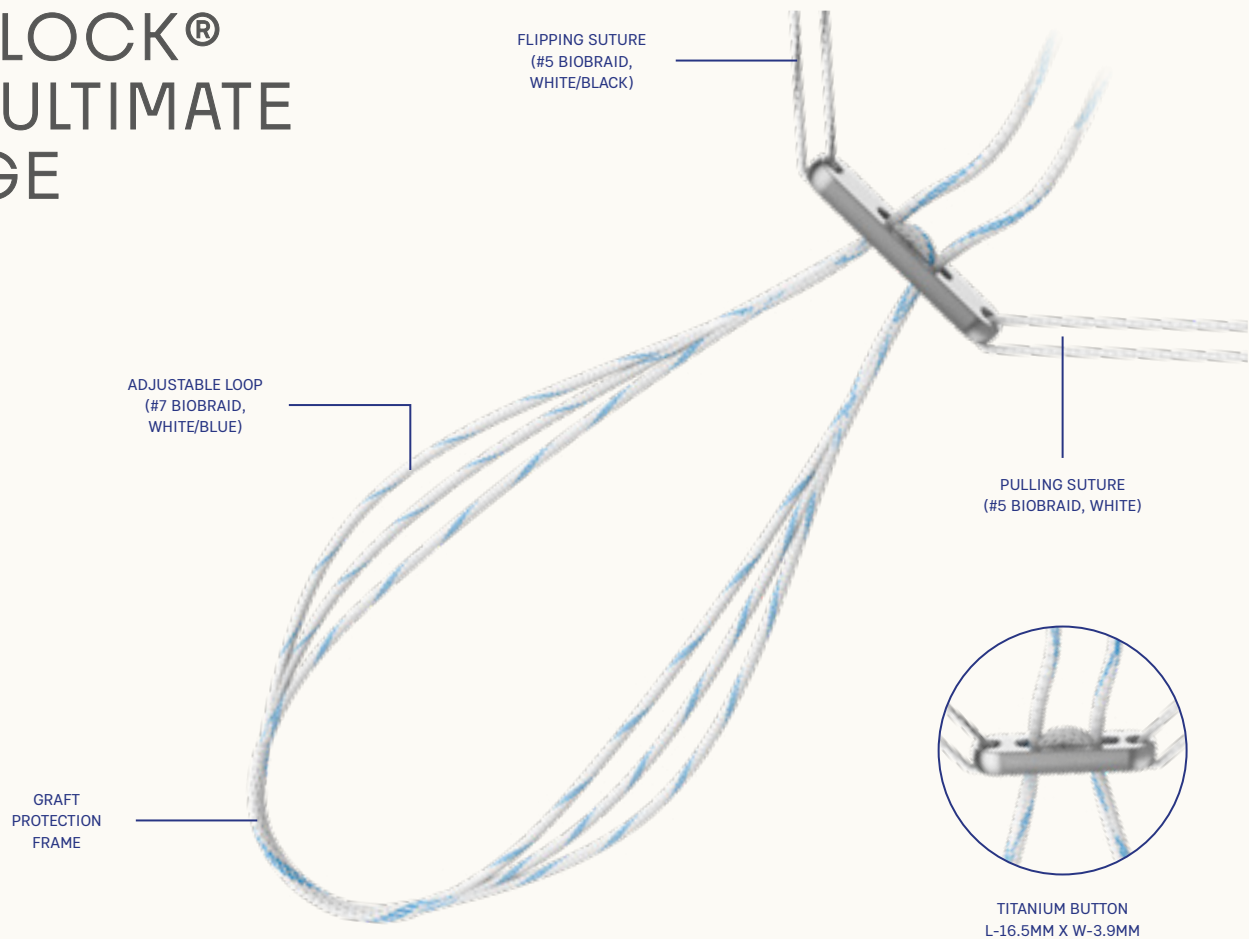

GFS Ultimate Large is an adjustable loop with an oblong shaped titanium button used for Cruciate Graft Reconstructions. GFS Ultimate is a knotless adjustable loop provides double-locking mechanism, which eliminates the need for knot tying. Because of the long length of the button GFS Ultimate Large is considered an excellent choice for revision surgeries.

GFS ULTIMATE LARGE PROVIDES THREE SUTURES (UHMWPE):

- **Adjustable Suture (#7, White/Blue):** It allows the surgeon to maximize the amount of graft inside the femoral tunnel, thereby optimizing the healing process. It also enables the calibration of the loop to its optimum size.
- **Pulling Suture (#5, White):** It is available to pull the graft inside the tunnel.
- **Flipping Suture (#5, White/Black):** It ensures the flipping of the button.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	LENGTH OF LOOP
6-006-02	AUXILOCK® GFS Ultimate Large Button, Button L: 16.5mm, W: 3.9mm (Adjustable Loop)	Package loop length: 70 mm Min: 12 mm

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-57	4.5mm Cannulated Headed Reamer