




AUXILOCK® GFS ULTIMATE MINI

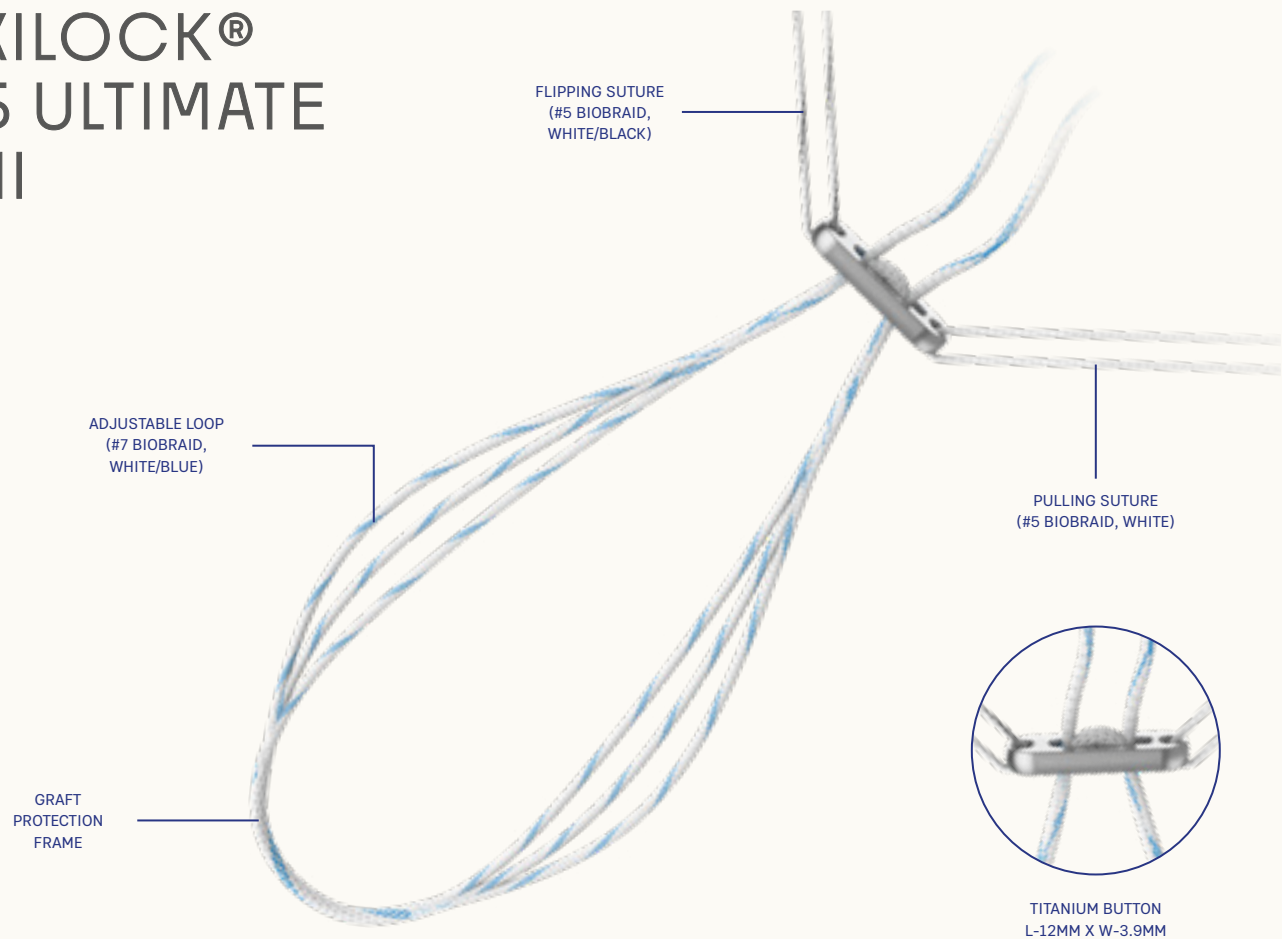

GFS Ultimate Mini is an adjustable loop with an oblong shaped titanium button used for Cruciate Graft Reconstructions. GFS Ultimate is a knotless adjustable loop provides double-locking mechanism, which eliminates the need for knot tying.

GFS ULTIMATE MINI PROVIDES THREE SUTURES (UHMWPE):

- **Adjustable Suture (#7, White/Blue):** It allows the surgeon to maximize the amount of graft inside the femoral tunnel, thereby optimizing the healing process. It also enables the calibration of the loop to its optimum size.
- **Pulling Suture (#5, White):** It is available to pull the graft inside the tunnel.
- **Flipping Suture (#5, White/Black):** It ensures the flipping of the button.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	LENGTH OF LOOP
6-006-01	AUXILOCK® GFS Ultimate Mini Button, Button L: 12mm, W: 3.9mm (Adjustable Loop)	Package loop length: 70 mm Min: 12 mm

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-57	4.5mm Cannulated Headed Reamer