




AUXILOCK® GFS ULTIMATE MINI BUTTON WITH BRACE

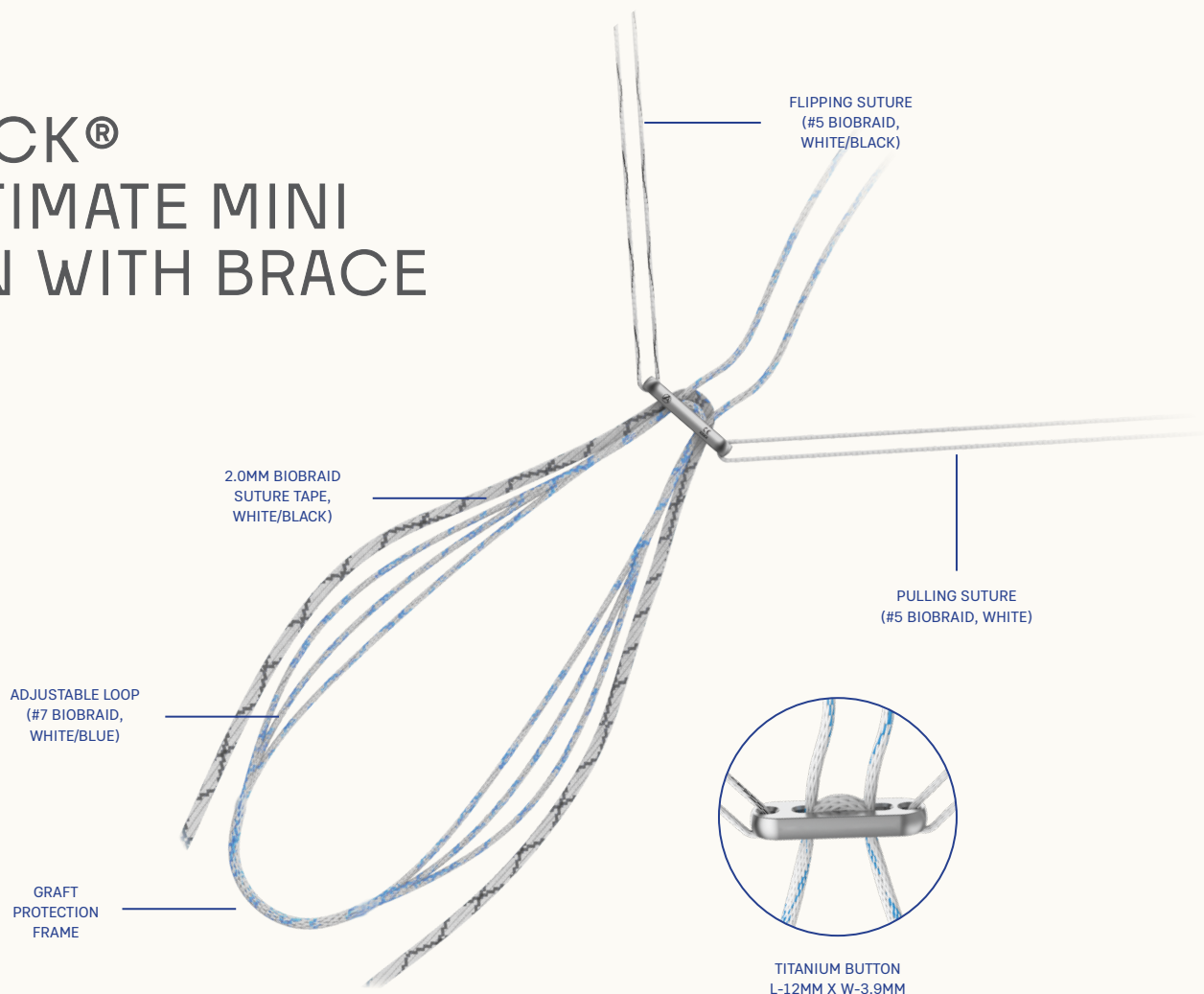

GFS ULTIMATE MINI BUTTON WITH BRACE implant delivers enhanced handling characteristics and superior resistance to graft abrasion and tissue pull-through compared to traditional round sutures. Designed with precision in mind, its adjustable-loop mechanism facilitates incremental re-tensioning of the graft construct even after cortical fixation, ensuring optimal tensioning and stability.

GFS Ultimate Mini provides four sutures (UHMWPE):

- **Adjustable Suture (#7, White/Blue):** It allows the surgeon to maximize the amount of graft inside the femoral tunnel, thereby optimizing the healing process. It also enables the calibration of the loop to its optimum size.
- **Pulling Suture (#5, White):** It is available to pull the graft inside the tunnel.
- **Flipping Suture (#5, White/Black):** It ensures the flipping of the button.
- **2.0mm BioBraid Suture Tape, White/black):** Offers more secured Tibial Fixation.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	LENGTH OF LOOP
6-006-05	AUXILOCK GFS Ultimate Mini Button With Brace, Button L: 12mm, W: 3.9mm (Adjustable Loop)	Package loop length: 70 mm Min: 12 mm
6-006-06	AUXILOCK® GFS Ultimate Large Button With Brace, Button L: 16.5mm, W: 3.9mm (Adjustable Loop)	Package loop length: 70 mm Min: 12 mm

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-57	4.5mm Cannulated Headed Reamer