




AUXILOCK® GYROLOCK PEEK OPTIMA KNOTLESS SCREW IN ANCHOR



SELF-PUNCHING
EYELET

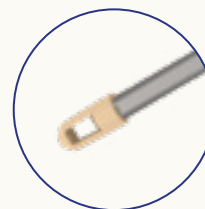

PEEK OPTIMA
EYELET



PEEK OPTIMA
ANCHOR

GyroLock Knotless PEEK OPTIMA Screw-in Anchors consists of PEEK-OPTIMA (Poly-ether-ether-ketone) material. The GyroLock suture anchor comes with Dual start thread technology (fast inserting suture anchor) which allows anchor deployment in 3.5 rotation.

Vented and cannulated GyroLock anchors are designed for repair and reconstruction with soft-tissue graft. Vented design promotes bone growth. It is also available in a self-punching eyelet, designed for the lateral row in RC repair.

The GyroLock is designed for multiple procedures in single or double row for lateral row applications. GyroLock suture anchor eyelet can accommodate up to four 1.4mm tapes or four #2 sutures.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SILICON COLOR
6-052-02	AUXILOCK® GyroLock 4.75mm x 22.5mm Knotless PEEK OPTIMA Screw-In Anchor	
6-052-03	AUXILOCK® Gyrolock 5.5mm x 22.5mm Knotless PEEK OPTIMA Screw-In Anchor	
6-052-04	AUXILOCK® GyroLock SP 4.75mm x 24.5mm Knotless PEEK OPTIMA Screw-In Anchor, Self-Punching	
6-052-05	AUXILOCK® GyroLock SP 5.5mm x 24.5mm Knotless PEEK OPTIMA Screw-In Anchor, Self-Punching	

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-69	Awl/Tap for 4.75/5.5mm GyroLock Screw-In Anchor
7-013-70	Tap for 4.75/5.5mm GyroLock Screw-In Anchor