





AUXILOCK® ROTADOR PEEK OPTIMA SCREW-IN ANCHOR

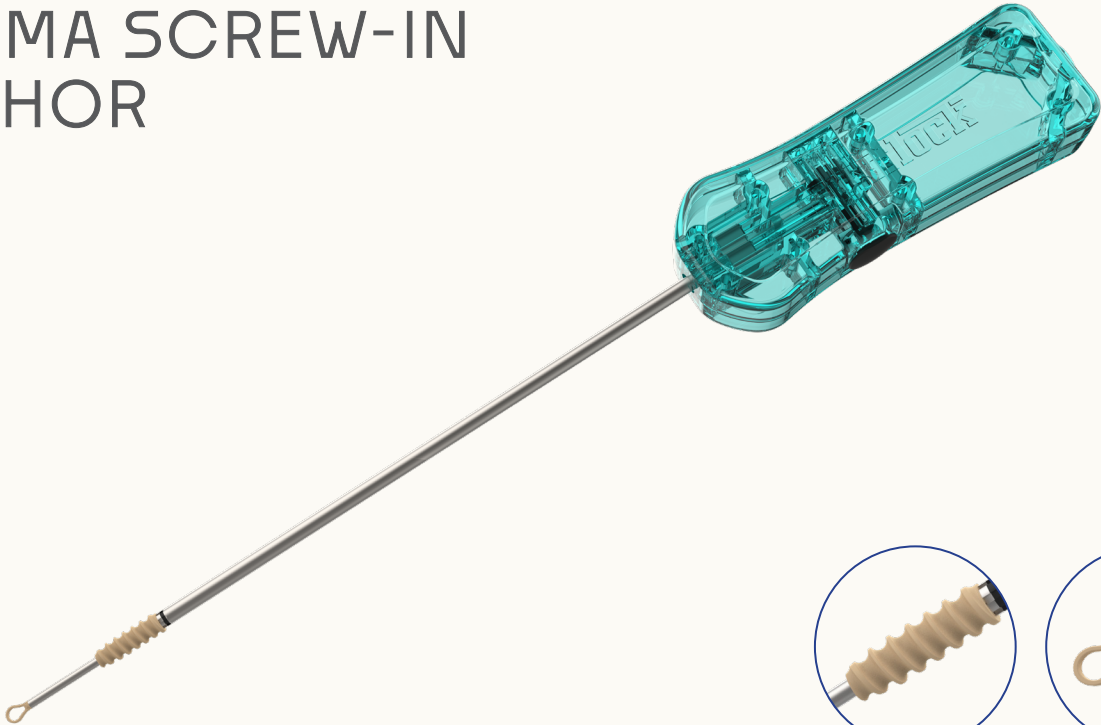

PEEK OPTIMA
ANCHOR

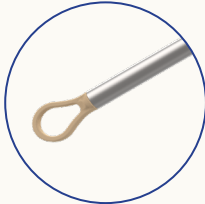

PEEK OPTIMA
EYELET

AUXILOCK® ROTADOR PEEK OPTIMA Screw-In Anchors are fully threaded knotless anchors available in 4.75, 5.5 and 6.25mm diameters with PEEK OPTIMA anchor body and eyelet. These anchors are designed to be used with sutures or tapes for rotator cuff repair using the ‘bridge’ technique.

Moreover, the ‘knotless’ technique consists of passing sutures or tapes of the medial row anchors through the tissue. They are finally inserted into the bone socket once they’re loaded through the Rotador anchor eyelet. This technique eliminates possible complications caused by knots compared to other conventional anchors. The eyelet of the anchor accommodates up to 6 sutures and 2 tapes at a time.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-020-01	AUXILOCK® Rotador 4.75mm x 15mm PEEK OPTIMA Screw-In Anchor
6-020-02	AUXILOCK® Rotador 5.5mm x 15mm PEEK OPTIMA Screw-In Anchor
6-020-03	AUXILOCK® Rotador 6.25mm x 15mm PEEK OPTIMA Screw-In Anchor

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-24	Ø4.1mm Awl	6-020-01, 6-020-02, 6-020-03