




AUXILOCK® SYNDESMOSIS ULTIMATE KIT

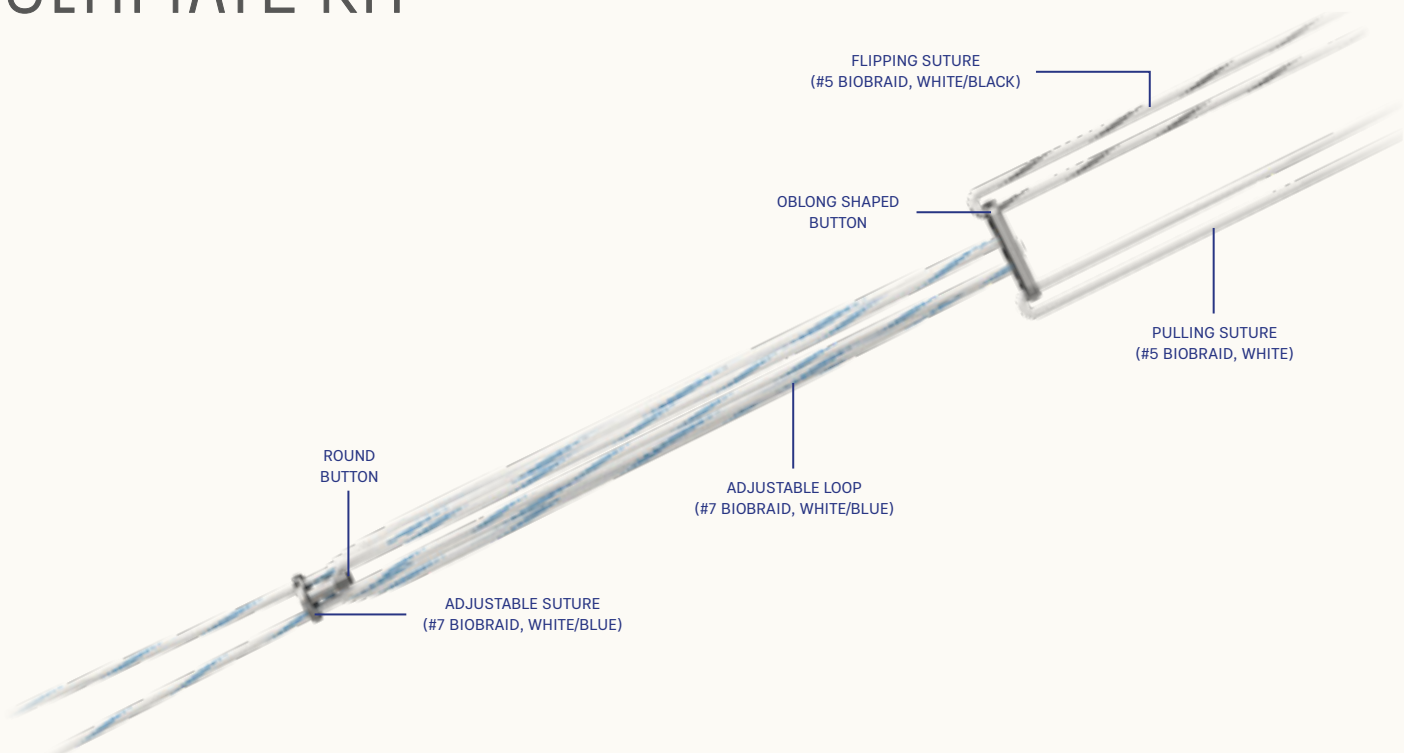

AUXILOCK Syndesmosis Ultimate Kit can be used for the fixation of Ankle ® Syndesmosis joint. The versatile design of the kit is ideal for use with or without a fracture plate. The kit comprises of adjustable loop with round shaped button (Dia. 5.5mm) and oblong shaped button (L 12mm X W 3.2mm X H 1.5mm).

The round button has a tapered slot design which can easily fix on any fracture plate. The kit is passed through 3.5mm drilling tunnel from lateral to medial provides the fixation of oblong shaped button on the medial side and round button on the lateral side of the syndesmosis joint. The syndesmosis kit provides a double locking mechanism which eliminates the need for knot tying.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-041-01	AUXILOCK® Syndesmosis Ultimate Kit

RELATED INSTRUMENTATION (INCLUDE IN KIT)

CODE	PRODUCT DESCRIPTION
7-147-01	Cannulated Drill Bit with Quick Coupling End, Ø3.5mm x Length 175mm
7-147-02	Ø1.5mm K-Wire
7-147-03	Ø2.4mm Passing Pin
7-147-04	Suture Pull Handle