




AUXILOCK® TENDO PEEK OPTIMA SCREW-IN ANCHOR



PEEK OPTIMA
WASHER



PEEK OPTIMA
TENDO SCREW

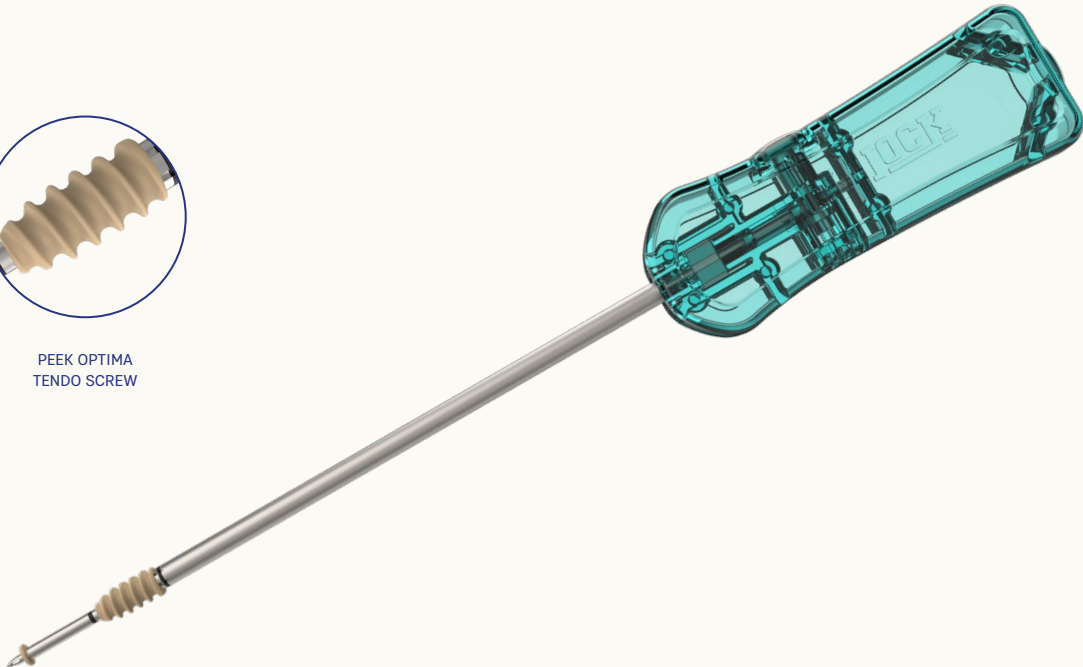

AUXILOCK® Tendo PEEK OPTIMA Screw-In Anchor is ideal for proximal and distal biceps tenodesis repair. It is available in three diameters 7.0, 8.0 & 9.0mm & two lengths 10mm and 15mm, with PEEK OPTIMA interference/tenodesis screws and a PEEK OPTIMA washer.

The Eyelet is specifically designed to hold the biceps tendon at the bottom of the drill hole. After fixing up the tendon into the drill hole, it can be locked into the bone by advancing the tenodesis screw. The anchor is designed to save steps and minimize the length of the procedure.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-042-01	AUXILOCK® Tendo 7.0mm x 15mm PEEK OPTIMA Screw-In Anchor
6-042-02	AUXILOCK® Tendo 8.0mm x 15mm PEEK OPTIMA Screw-In Anchor
6-042-03	AUXILOCK® Tendo 9.0mm x 15mm PEEK OPTIMA Screw-In Anchor
6-042-04	AUXILOCK® Tendo 7.0mm x 10mm PEEK OPTIMA Screw-In Anchor
6-042-05	AUXILOCK® Tendo 8.0mm x 10mm PEEK OPTIMA Screw-In Anchor
6-042-06	AUXILOCK® Tendo 9.0mm x 10mm PEEK OPTIMA Screw-In Anchor

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-58	Graft Sizer, 6.0mm to 9.0mm (1.0mm Steps)
7-013-59	Ø6.0mm Trocar Tip Headed Reamer
7-013-60	Ø7.0mm Trocar Tip Headed Reamer
7-013-61	Ø8.0mm Trocar Tip Headed Reamer
7-013-62	Ø9.0mm Trocar Tip Headed Reamer