




AUXILOCK® KNOTLESS PEEK CF PUSH-IN SUTURE ANCHOR

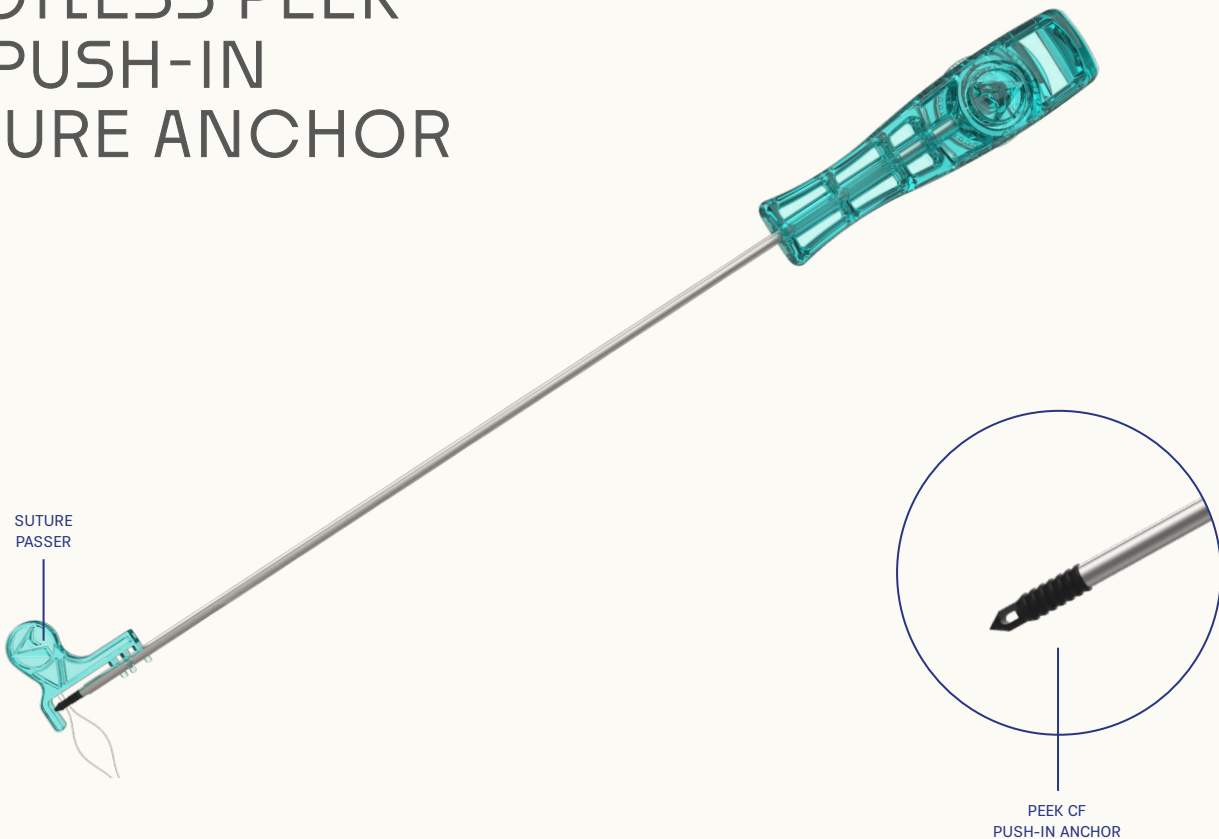

AUXILOCK® knotless push-in anchors are a step-saving alternative to conventional knotted suture anchors. AUXILOCK® knotless push-in anchors are recommended for use in Bankart and SLAP repair surgeries.

These Anchors are made of carbon fiber-reinforced PEEK-OPTIMA and are available in a variety of diameters of 2.8mm and 3.5mm providing intraoperative flexibility. The PEEK CF provides modulus of elasticity closely matching the cortical bone. It is non-absorbable, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-003-01	AUXILOCK® 2.8mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer
6-003-02	AUXILOCK® 3.5mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-003-01, 6-003-02
7-013-54	Ø2.5mm X 200 mm Drill Bit (For Knotless Anchor)	6-003-01
7-013-55	Ø3.2mm X 200 mm Drill Bit (For Knotless Anchor)	6-003-02
7-013-02	8 Point Drill Guide, Length 100mm	6-003-01, 6-003-02

SHOULDER INSTABILITY REPAIR