





SHOULDER



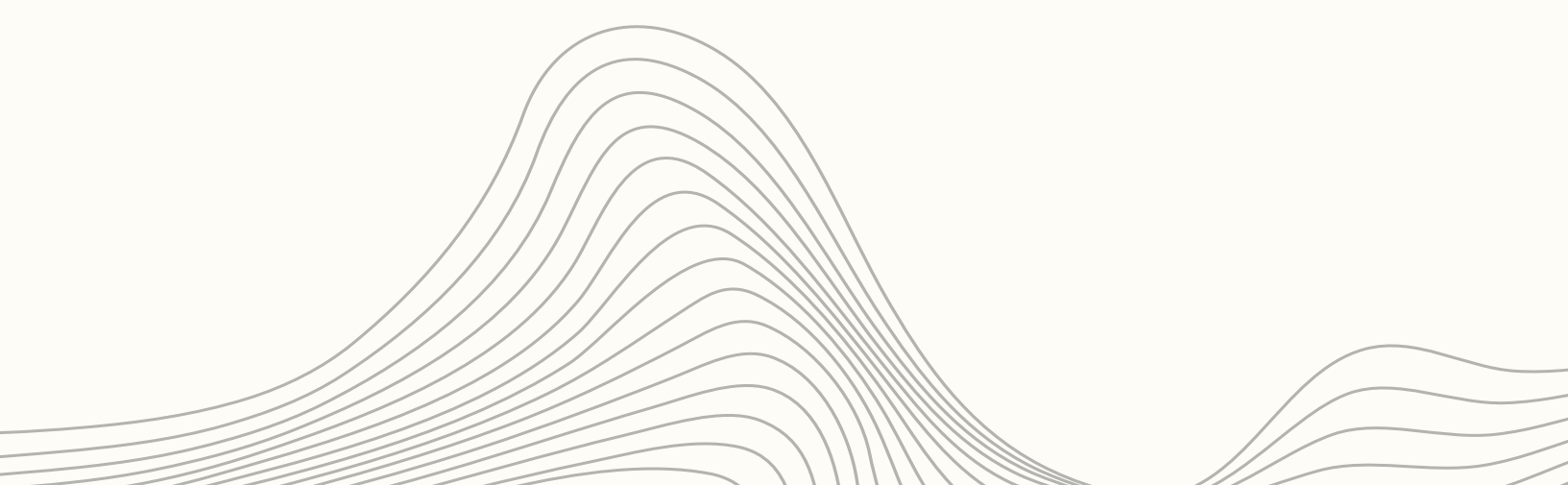

AUXILOCK®AUXPIK SUTURE-BASED ANCHOR

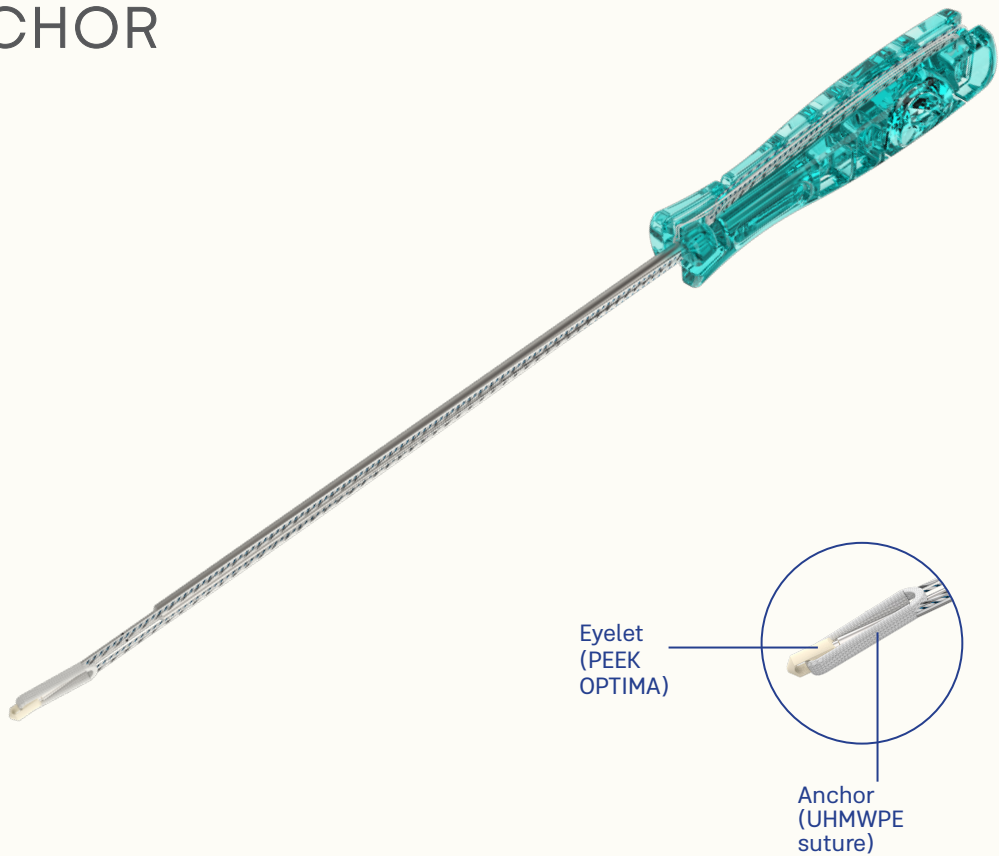

AUXILOCK® AuxPik suture-based anchor is made with UHMWPE suture anchor body and PEEK OPTIMA eyelet tip.

The anchor is designed to deliver efficiency and promote ease of use. The AuxPik suture-based anchor is available in 3.2mm diameter with various combination of BioBraid and suture tapes. It can be used for rotator cuff repair surgeries. The AUXILOCK® AuxPik suture-based anchor provides a small footprint and also asserts subcortical fixation for anchor insertion. The drill guide and anchor driver combination is well designed to improve the performance and reliability.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR
6-001-03	AUXILOCK® 3.2mm AuxPik Suture-Based Anchor With Two #2 Bio-Braid: White/Blue & White/Black	
6-001-09	AUXILOCK® 3.2mm AuxPik Suture-Based Anchor With One #2 BioBraid White/Black and One 1.4mm Suture Tape:White/Blue	

RECOMMENDED INSTRUMENTATION


CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-001-03, 6-001-09
7-013-13	Ø3.2mm X 200mm Drill Bit	6-001-03, 6-001-09
7-013-23	Ø3.2mm Awl	6-001-03, 6-001-09
7-013-02	8 Point Drill Guide, Length 100mm	6-001-03, 6-001-09

ROTATOR CUFF REPAIR



AUXILOCK® AUXSUTURE ANCHOR

ROTATOR CUFF REPAIR

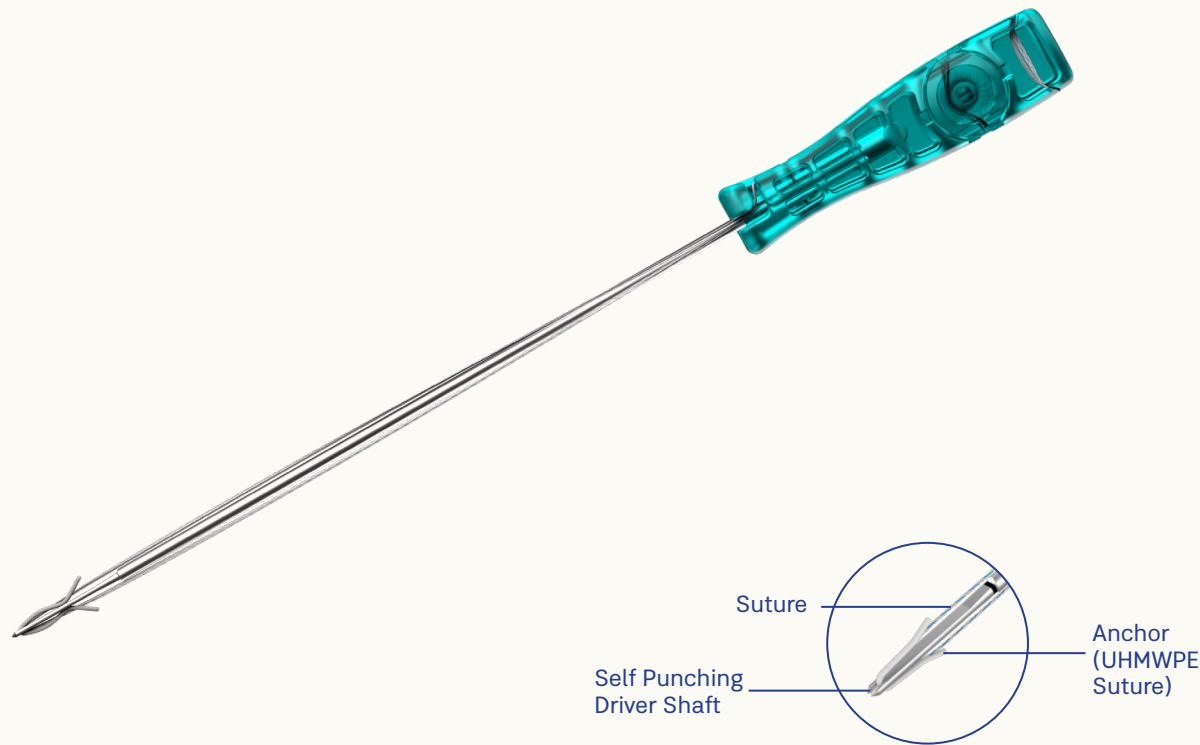

AUXILOCK® AuxSuture Anchor provides numerous advantages over traditional anchors. The strong fixation, less bone removal and a push in technique provide higher efficiency in rotator cuff repair procedures.

AuxSuture anchor is available in diameter of 2.9mm. It is a self-punching anchor that can be inserted without a pilot hole in the cancellous bone. The specially designed 2.9mm Awl for AuxSuture Anchor provides an accurate hole placement for the anchor in the hard bones. The small 2.9mm footprint of AuxSuture anchor provides placement and fixation advantages when bone real estate is limited, especially during revision surgeries. AuxSuture Anchor is available in double and triple loaded options with various combination of BioBraid and suture tapes.

RECOMMENDED INSTRUMENTATION


CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-22	Ø2.9mm Awl	6-036-01, 6-036-02, 6-036-03, 6-036-04, 6-036-05

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-036-01	AUXILOCK® 2.9mm AuxSuture Anchor With Two #2 BioBraid: White/Blue & White/Green		
6-036-02	AUXILOCK® 2.9mm AuxSuture Anchor With Two #2 BioBraid: White/Blue & White/Green,With Needle:M0-6		
6-036-03	AUXILOCK® 2.9mm AuxSuture Anchor With Three #2 BioBraid: White, White/Blue & White/Green		
6-036-04	AUXILOCK® 2.9mm AuxSuture Anchor With One #2 BioBraid White/Black and One 1.4mm Suture Tape: White/Blue		
6-036-05	AUXILOCK® 2.9mm AuxSuture Anchor With Two 1.4mm Suture Tape: White/Black & White/Blue		



AUXILOCK® PEEK CF SCREW-IN SUTURE ANCHOR

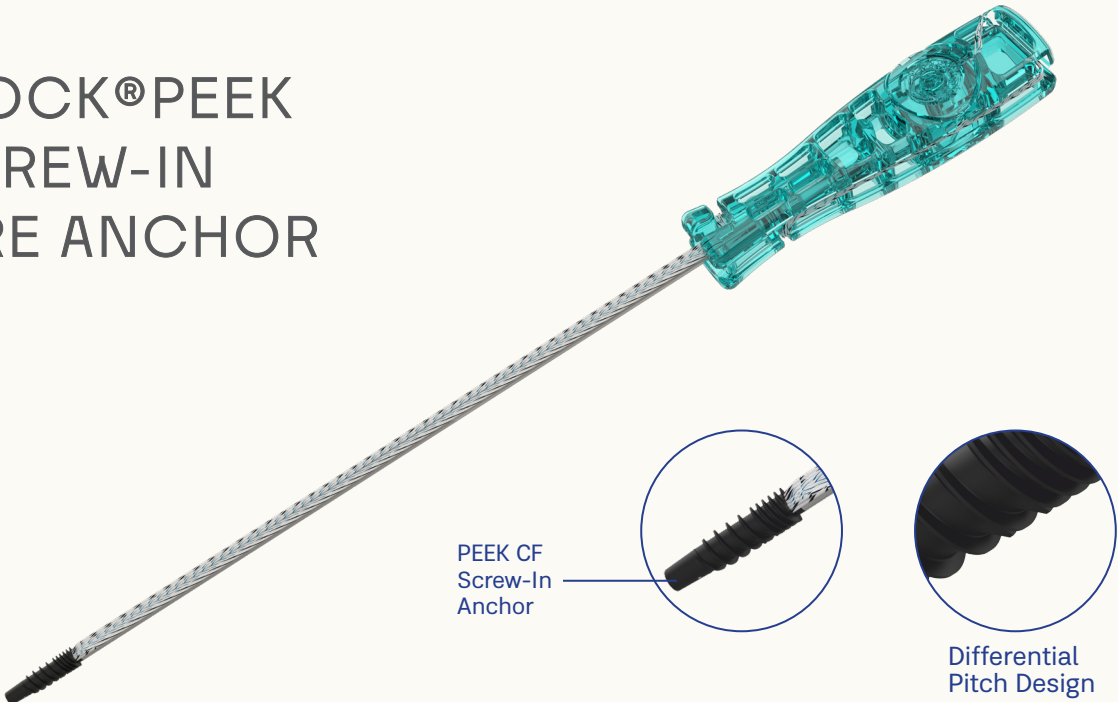


ROTATOR CUFF REPAIR

AUXILOCK® PEEK CF Screw-In Suture anchor is a fully threaded suture anchor featuring dual threads to maximize cortical and cancellous fixation. Differential pitch design improves fixation by optimizing stress distribution between dense cortical and weaker cancellous bone. This design minimizes stress concentrations, reduces the risk of anchor pullout, and enhances stability in patients with poor bone quality by maximizing surface area engagement. It is available in diameter of 4.5, 5.5 and 6.5mm with double and triple loaded options. It has a flat tip to facilitate the insertion. The anchor is particularly suitable for repairing rotator cuff and associated pathologies

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR
6-046-01	AUXILOCK® 4.5mm PEEK CF Screw-In Suture Anchor With Three #2 BioBraid: White, White/Blue & White/Black	
6-046-02	AUXILOCK® 4.5mm PEEK CF Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black	
6-046-03	AUXILOCK® 5.5mm PEEK CF Screw-In Suture Anchor With Three #2 BioBraid: White, White/Blue & White/Black	

6-046-04	AUXILOCK® 5.5mm PEEK CF Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black	
6-046-05	AUXILOCK® 6.5mm PEEK CF Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black	
6-046-06	AUXILOCK® 4.5mm PEEK CF Screw-In Suture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Black & White/Blue	
6-046-07	AUXILOCK® 5.5mm PEEK CF Screw-In Suture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Black & White/Blue	
6-046-08	AUXILOCK® 6.5mm PEEK CF Screw-In Suture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Black & White/Blue	
6-046-09	AUXILOCK® 4.5mm PEEK CF Screw-In Suture Anchor With Two 1.4mm Suture Tape: White/Blue & White/Black	
6-046-10	AUXILOCK® 5.5mm PEEK CF Screw-In Suture Anchor With Two 1.4mm Suture Tape: White/Blue & White/Black	
6-046-11	AUXILOCK® 6.5mm PEEK CF Screw-In Suture Anchor With Two 1.4mm Suture Tape: White/Blue & White/Black	

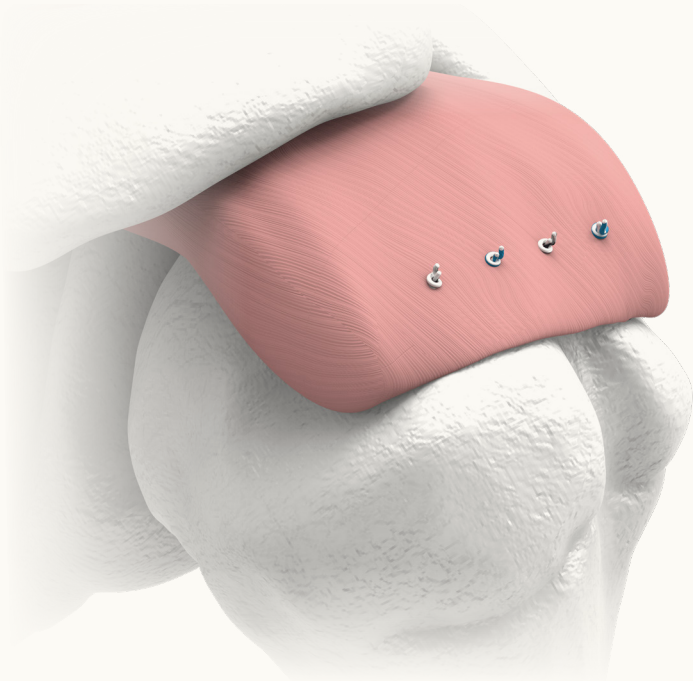

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-016-21	AUXILOCK® 4.5mm PEEK CF Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-016-22	AUXILOCK® 5.5mm PEEK CF Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-016-23	AUXILOCK® 6.5mm PEEK CF Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	BONE DENSITY	COMPATIBLE IMPLANTS
7-013-25	Awl For 4.5/5.5mm PEEK Screw-In Anchor	Soft Bone	6-046-01, 6-046-02, 6-046-03, 6-046-04, 6-046-06, 6-046-07, 6-046-09, 6-046-10, 6-046-21, 6-046-22
7-013-26	Awl For 6.5mm PEEK Screw-In Anchor		6-046-05, 6-046-08, 6-046-11, 6-046-23
7-013-27	Awl/ Tap For 4.5/5.5mm PEEK Screw-In Anchor	Normal Bone	6-046-01, 6-046-02, 6-046-03, 6-046-04, 6-046-06, 6-046-07, 6-046-09, 6-046-10, 6-046-21, 6-046-22
7-013-28	Awl/ Tap For 6.5mm PEEK Screw-In Anchor		6-046-05, 6-046-08, 6-046-11, 6-046-23
7-013-29	Tap For 4.5/5.5mm PEEK Screw-In Anchor	Hard Bone	6-046-01, 6-046-02, 6-046-03, 6-046-04, 6-046-06, 6-046-07, 6-046-09, 6-046-10, 6-046-21, 6-046-22
7-013-30	Tap For 6.5mm PEEK Screw-In Anchor		6-046-05, 6-046-08, 6-046-11, 6-046-23



AUXILOCK® PEEK OPTIMA SCREW-IN SUTURE ANCHOR

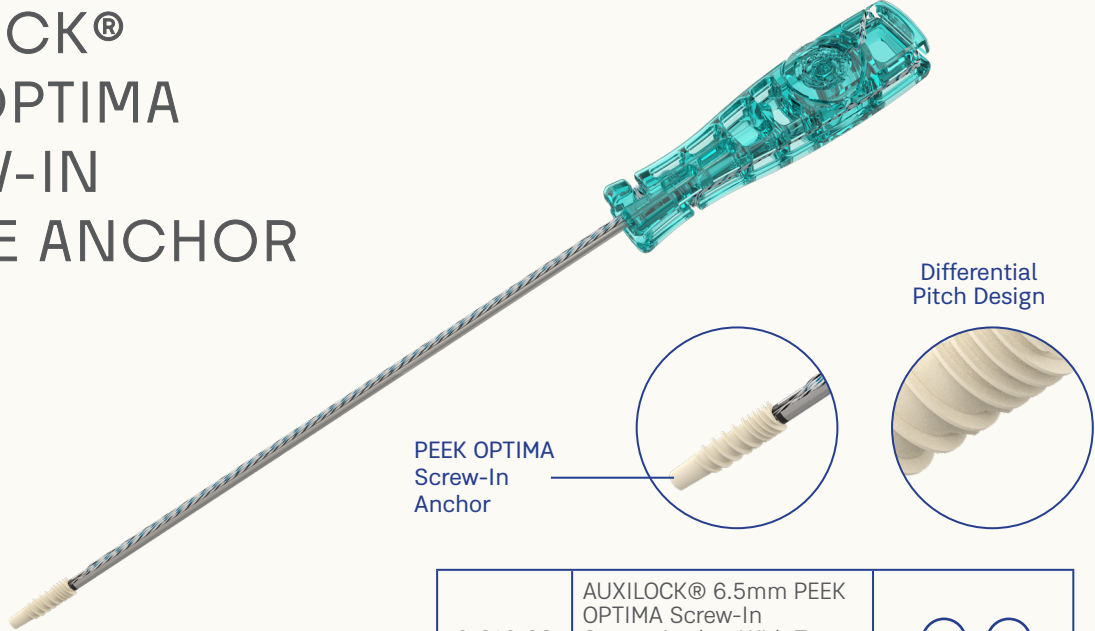


AUXILOCK® PEEK OPTIMA Screw-In Suture anchor is a fully threaded suture anchor featuring dual threads to maximize cortical and cancellous fixation. Differential pitch design improves fixation by optimizing stress distribution between dense cortical and weaker cancellous bone. This design minimizes stress concentrations, reduces the risk of anchor pullout, and enhances stability in patients with poor bone quality by maximizing surface area engagement. It is available in diameter of 4.5, 5.5 and 6.5mm with double and triple loaded options. It has a flat tip to facilitate the insertion. The anchor is particularly suitable for repairing rotator cuff and associated pathologies

It has a flat tip to facilitate the insertion. The anchor is particularly suitable for repairing rotator cuff and associated pathologies. The anchors are also available with needles which are ideal for mini-open rotator cuff repair procedures.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR
6-016-05	AUXILOCK® 4.5mm PEEK OPTIMA Screw-In Suture Anchor With Two #2 BioBraid: White/ Blue & White/Black	
6-016-07	AUXILOCK® 5.5mm PEEK OPTIMA Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black	

6-016-08	AUXILOCK® 6.5mm PEEK OPTIMA Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black	
6-016-13	AUXILOCK® 4.5mm PEEK OPTIMA Screw-In Suture Anchor With Three #2 BioBraid: White, White/Blue & White/Black	
6-016-14	AUXILOCK® 5.5mm PEEK OPTIMA Screw-In Suture Anchor With Three #2 BioBraid: White,White/ Blue & White/Black	
6-016-15	AUXILOCK® 4.5mm PEEK OPTIMA Screw-In Suture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Blue & White/Black	
6-016-16	AUXILOCK® 5.5mm PEEK OPTIMA Screw-In Suture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Blue & White/Black	
6-016-17	AUXILOCK® 6.5mm PEEK OPTIMA Screw-In Suture Anchor With One #2 Bio-Braid and One 1.4mm Suture Tape: White/Blue & White/Black	
6-016-18	AUXILOCK® 4.5mm PEEK OPTIMA Screw-In Suture Anchor With Two 1.4mm Suture Tape: White/Blue & White/Black	
6-016-19	AUXILOCK® 5.5mm PEEK OPTIMA Screw-In Suture Anchor With Two 1.4mm Suture Tape: White/Blue & White/Black	
6-016-20	AUXILOCK® 6.5mm PEEK OPTIMA Screw-In Suture Anchor With Two 1.4mm Suture Tape: White/Blue & White/Black	

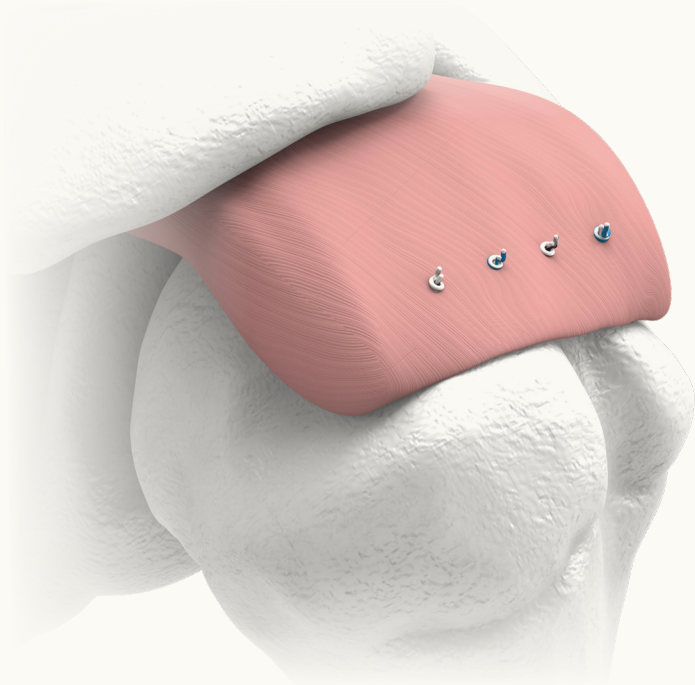

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-016-01	AUXILOCK® 4.5mm PEEK OPTIMA Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-016-03	AUXILOCK® 5.5mm PEEK OPTIMA Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-016-04	AUXILOCK® 6.5mm PEEK OPTIMA Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	BONE DENSITY	COMPATIBLE IMPLANTS
7-013-25	Awl For 4.5/5.5mm PEEK Screw-In Anchor	Soft Bone	6-016-05, 6-016-07, 6-016-13, 6-016-14, 6-016-15, 6-016-16, 6-016-18, 6-016-19, 6-016-01, 6-016-03
7-013-26	Awl For 6.5mm PEEK Screw-In Anchor		6-016-08, 6-016-17, 6-016-20, 6-016-04
7-013-27	Awl/ Tap For 4.5/5.5mm PEEK Screw-In Anchor	Normal Bone	6-016-05, 6-016-07, 6-016-13, 6-016-14, 6-016-15, 6-016-16, 6-016-18, 6-016-19, 6-016-01, 6-016-03
7-013-28	Awl/ Tap For 6.5mm PEEK Screw-In Anchor		6-016-08, 6-016-17, 6-016-20, 6-016-04
7-013-29	Tap For 4.5/5.5mm PEEK Screw-In Anchor	Hard Bone	6-016-05, 6-016-07, 6-016-13, 6-016-14, 6-016-15, 6-016-16, 6-016-18, 6-016-19, 6-016-01, 6-016-03
7-013-30	Tap For 6.5mm PEEK Screw-In Anchor		6-016-08, 6-016-17, 6-016-20, 6-016-04



AUXILOCK® TITANIUM SCREW-IN SUTURE ANCHOR

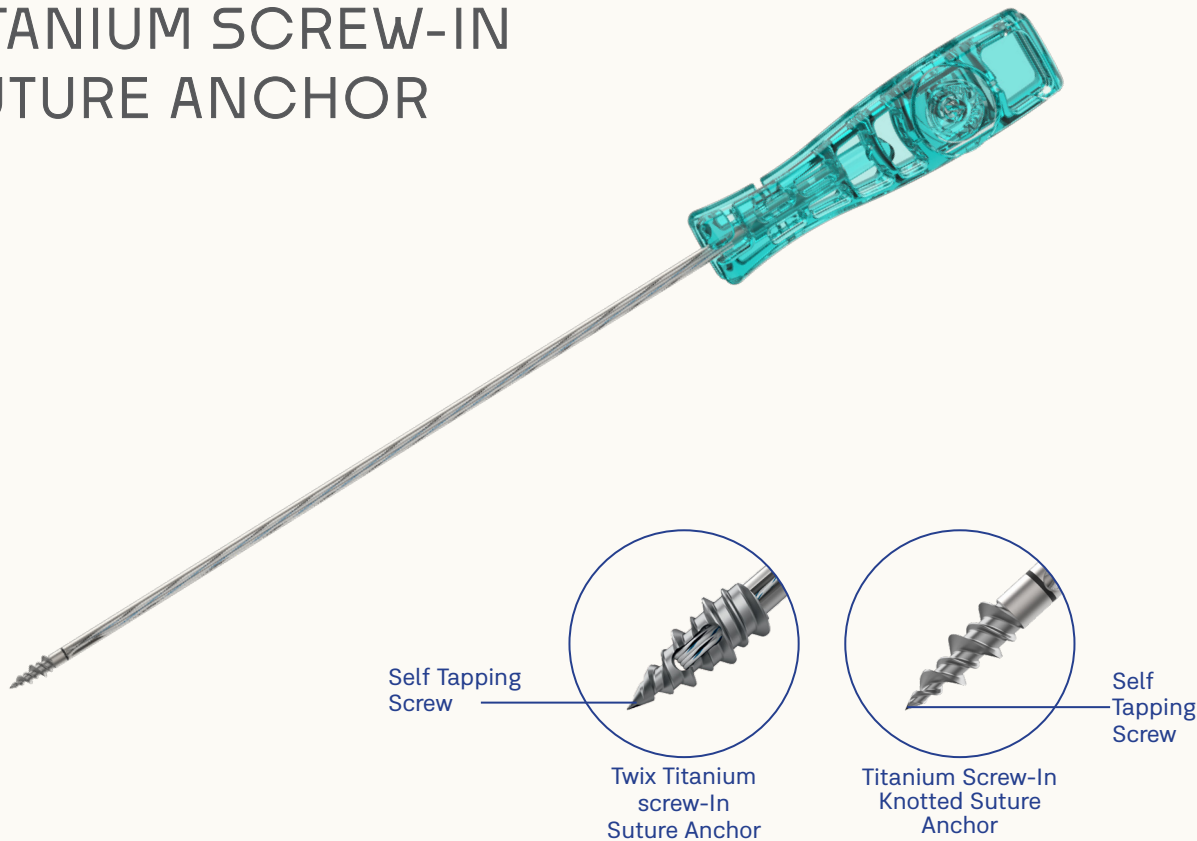

The self-tapping AUXILOCK® Titanium screw-in suture anchor is mostly used as a medial row anchor in rotator cuff repair surgery. It is recommended for use in small and large-joint repairs without the need for tapping or drilling.

AUXILOCK® Titanium screw-in anchors are available in diameter of 5.0 and 6.5mm with double and triple loaded options. The anchor design incorporates a cancellous thread with a very small core diameter to maximize pull-out strength in cancellous or osteoporotic bone. The anchor is designed for ultimate mechanical properties (pull-out strength, tensile strength, etc.) and ease-of-use. The self tap & self screw anchors are also available with needles which are ideal for mini-open rotator cuff repair procedures.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-004-02	AUXILOCK® 5.0mm Titanium Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black		

6-004-13	AUXILOCK® 5.0mm Titanium Screw-In Suture Anchor With Three #2 BioBraid: White, White/Blue & White/Black		
6-004-04	AUXILOCK® 6.5mm Titanium Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black		
6-004-05	AUXILOCK® 5.0mm Titanium Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-004-06	AUXILOCK® 6.5mm Titanium Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-045-01	AUXILOCK® 5.0mm Twix Titanium Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-045-02	AUXILOCK® 6.5mm Twix Titanium Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		



AUXILOCK® KNOTLESS PEEK CF PUSH-IN SUTURE ANCHOR

ROTATOR CUFF REPAIR

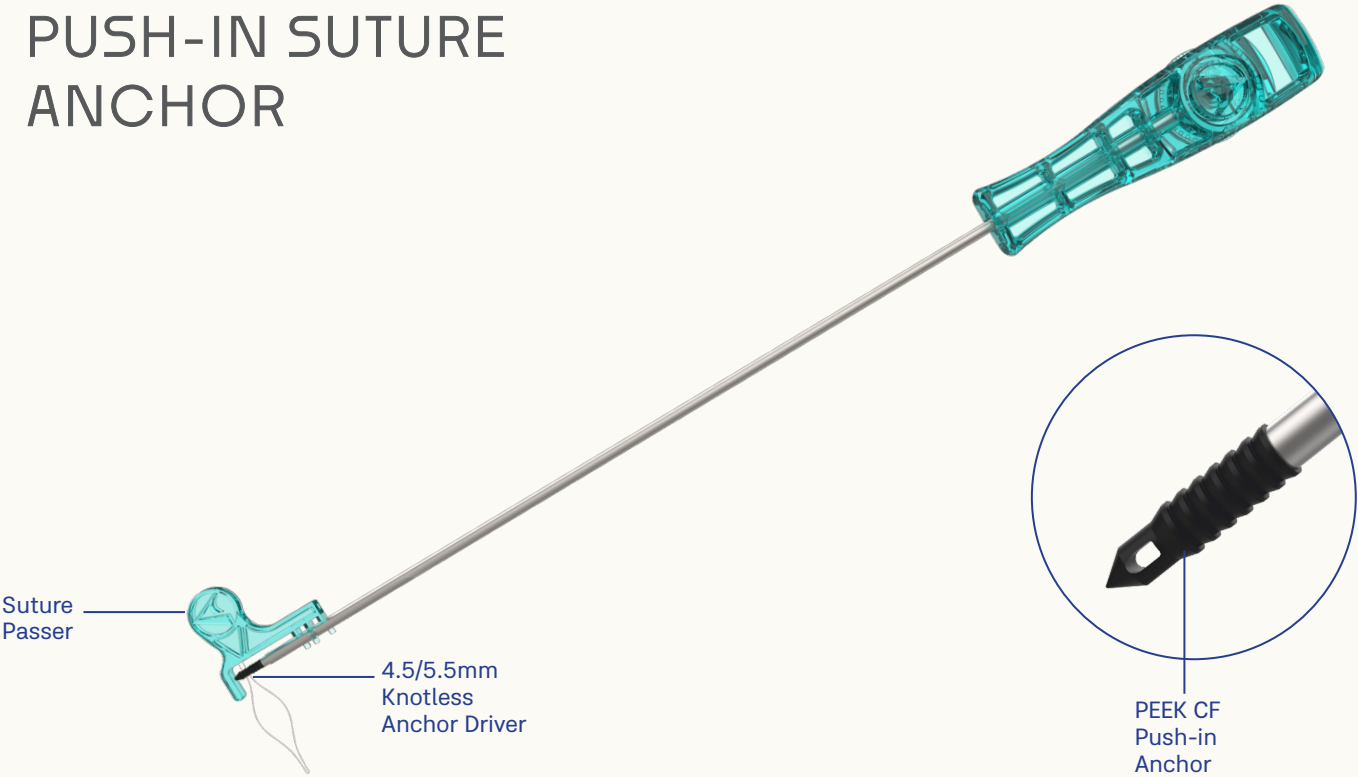

AUXILOCK® knotless push-in anchors provide a step-saving alternative to conventional knotted suture anchors. AUXILOCK® knotless push-in anchors are recommended for use in both large and small-joint soft tissue repairs. PEEK CF Push-In Suture Anchors consists of carbon fiber reinforced PEEK and are available in two diameters of 4.5mm and 5.5mm providing intraoperative flexibility.

PEEK CF provides modulus of elasticity closely matching the cortical bone. It is non-absorbable, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation. The anchor is available with and without driver. The 4.5/5.5mm knotless anchor driver is provided separately in the instrument set. It is reusable driver that minimizes the cost of the surgery unlike the pre-loaded suture anchors.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION
6-003-04	AUXILOCK® 4.5mm Knotless PEEK CF Push-In Suture Anchor With Suture Passer
6-003-05	AUXILOCK® 5.5mm Knotless PEEK CF Push-In Suture Anchor With Suture Passer
6-003-14	AUXILOCK® 4.5mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer
6-003-15	AUXILOCK® 5.5mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-23	Ø3.2mm Awl	6-003-04, 6-003-14
7-013-24	Ø4.1mm Awl	6-003-05, 6-003-15
7-013-32	4.5/5.5mm Knotless PEEK Anchor Driver	6-003-04, 6-003-05



AUXILOCK® KNOTLESS PEEK OPTIMA PUSH-IN SUTURE ANCHOR

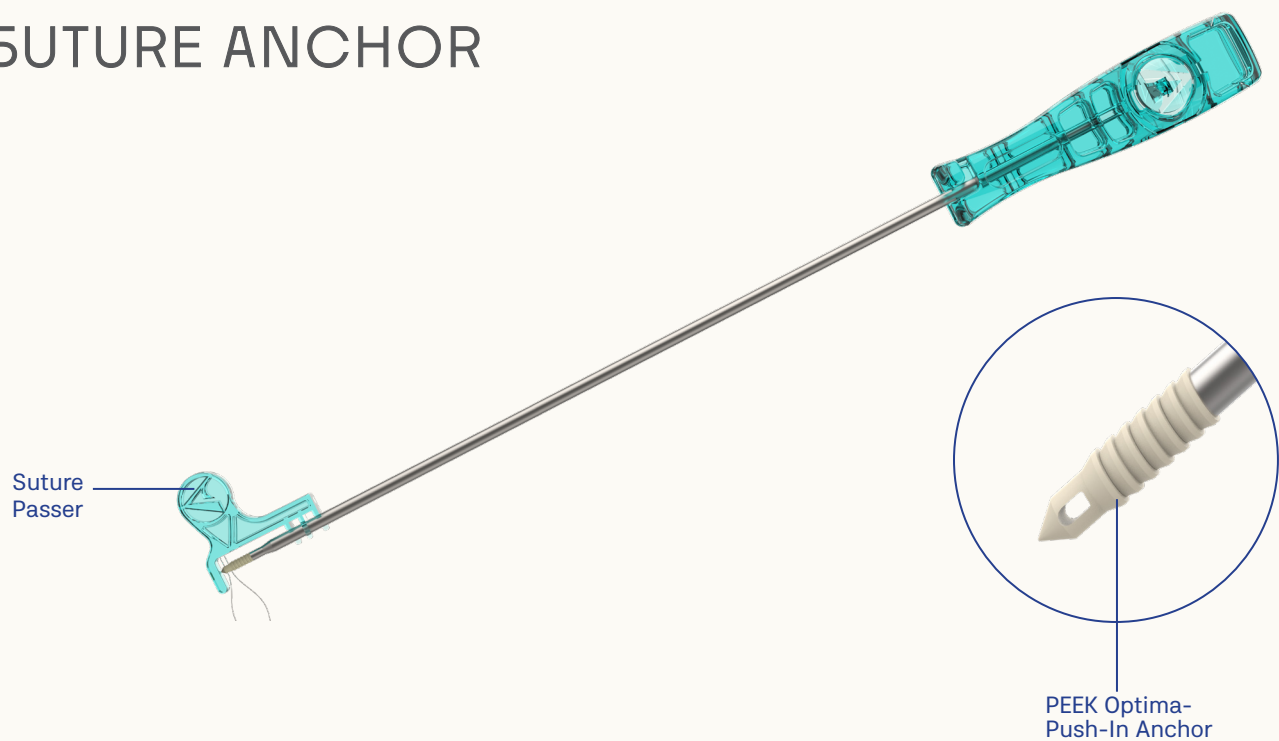

ROTATOR CUFF REPAIR

AUXILOCK® knotless Push-in anchors provide a step-saving alternative to conventional knotted suture anchors. AUXILOCK® knotless push-in anchors are recommended for use in both large and small-joint soft tissue repairs. PEEK OPTIMA Push-In Suture Anchors consists of PEEK OPTIMA (Poly-ether-ether-ketone) and are available in two diameters of 4.5mm and 5.5mm providing intraoperative flexibility.

The PEEK OPTIMA Push-In Suture Anchor is a non-absorbable PEEK suture anchor with an eyelet. Hence, providing superior abrasion resistance due to PEEK’s low coefficient friction. It is biocompatible, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation. The anchor is available with and without driver. The 4.5/5.5mm knotless anchor driver is provided separately in the instrument set. It is reusable driver that minimizes the cost of the surgery unlike the pre-loaded suture anchors.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION
6-031-04	AUXILOCK® 4.5mm Knotless PEEK OPTIMA Push-In Suture Anchor With Suture Passer
6-031-05	AUXILOCK® 5.5mm Knotless PEEK OPTIMA Push-In Suture Anchor With Suture Passer
6-031-07	AUXILOCK® 4.5mm Knotless PEEK OPTIMA Push-In Suture Anchor With Driver, With Suture Passer
6-031-08	AUXILOCK® 5.5mm Knotless PEEK OPTIMA Push-In Suture Anchor With Driver, With Suture Passer

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-23	Ø3.2mm Awl	6-031-04, 6-031-07
7-013-24	Ø4.1mm Awl	6-031-05, 6-031-08
7-013-32	4.5/5.5mm Knotless PEEK Anchor Driver	6-031-04, 6-031-05



AUXILOCK® ROTADOR PEEK CF SCREW-IN ANCHOR

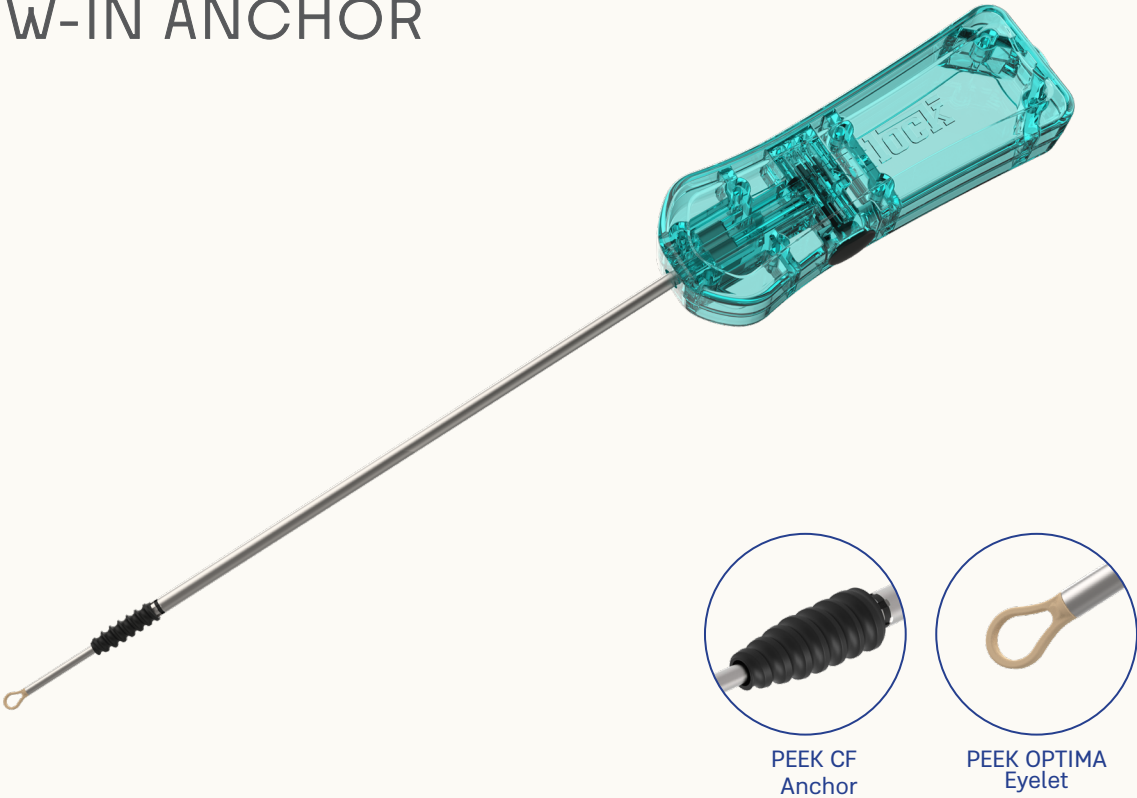

AUXILOCK® ROTADOR PEEK CF Screw-In Anchors are fully threaded knotless anchors available in 4.75, 5.5 and 6.25mm diameters with PEEK CF anchor body and PEEK OPTIMA eyelet.

These anchors are designed to be used with sutures or tapes for rotator cuff repair using the ‘bridge’ techniques. Moreover, the ‘knotless’ technique consists of passing sutures or tapes of the medial row anchors through the tissue. They are finally inserted into the bone socket once they’re loaded through the Rotador anchor eyelet. This technique eliminates possible complications caused by knots compared to other conventional anchors. The eyelet of the anchor accommodates up to 6 sutures and 2 tapes at a time.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION
6-020-04	AUXILOCK® Rotador 4.75mm x 15mm PEEK CF Screw-In Anchor
6-020-05	AUXILOCK® Rotador 5.5mm x 15mm PEEK CF Screw-In Anchor
6-020-06	AUXILOCK® Rotador 6.25mm x 15mm PEEK CF Screw-In Anchor

RECOMMENDED INSTRUMENTATION


CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-24	Ø4.1mm Awl	6-020-04, 6-020-05, 6-020-06



AUXILOCK® ROTADOR PEEK OPTIMA SCREW-IN ANCHOR

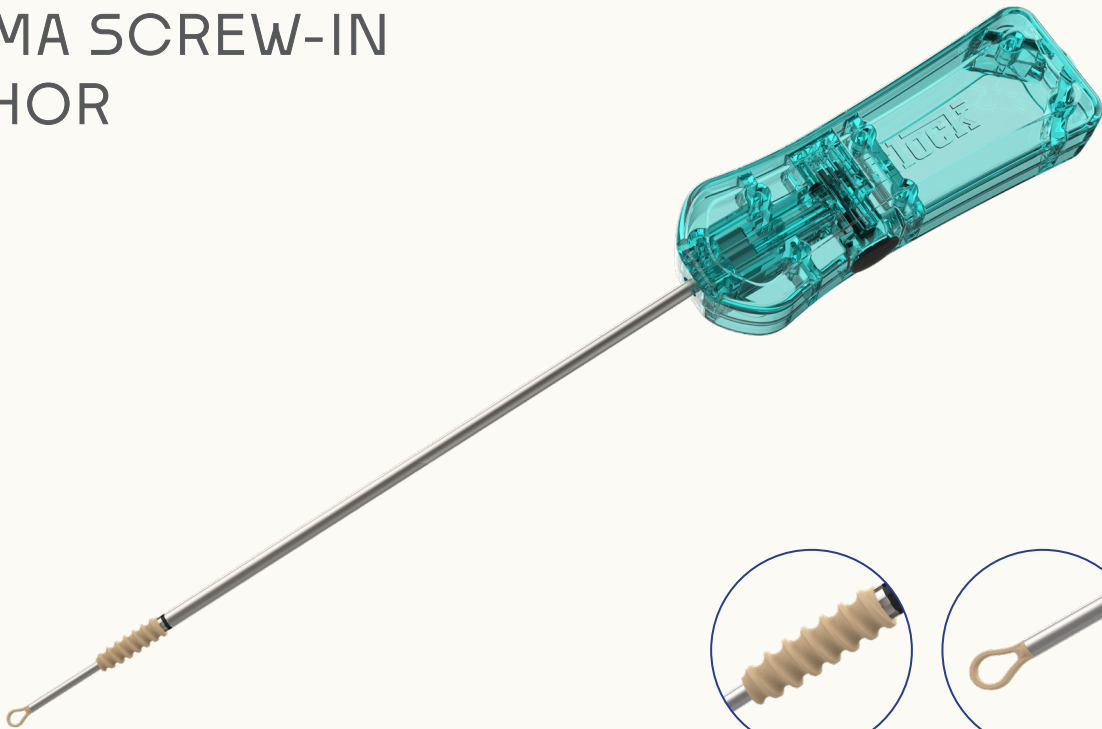

PEEK OPTIMA
Anchor

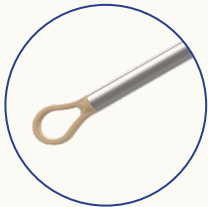

PEEK OPTIMA
Eyelet

AUXILOCK® ROTADOR PEEK OPTIMA Screw-In Anchors are fully threaded knotless anchors available in 4.75, 5.5 and 6.25mm diameters with PEEK OPTIMA anchor body and eyelet. These anchors are designed to be used with sutures or tapes for rotator cuff repair using the ‘bridge’ technique.

Moreover, the ‘knotless’ technique consists of passing sutures or tapes of the medial row anchors through the tissue. They are finally inserted into the bone socket once they’re loaded through the Rotador anchor eyelet. This technique eliminates possible complications caused by knots compared to other conventional anchors. The eyelet of the anchor accommodates up to 6 sutures and 2 tapes at a time.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-020-01	AUXILOCK® Rotador 4.75mm x 15mm PEEK OPTIMA Screw-In Anchor
6-020-02	AUXILOCK® Rotador 5.5mm x 15mm PEEK OPTIMA Screw-In Anchor
6-020-03	AUXILOCK® Rotador 6.25mm x 15mm PEEK OPTIMA Screw-In Anchor


RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-24	Ø4.1mm Awl	6-020-01, 6-020-02, 6-020-03


ROTATOR CUFF REPAIR



AUXILOCK® GYROLOCK PEEK OPTIMA KNOTLESS SCREW-IN ANCHOR

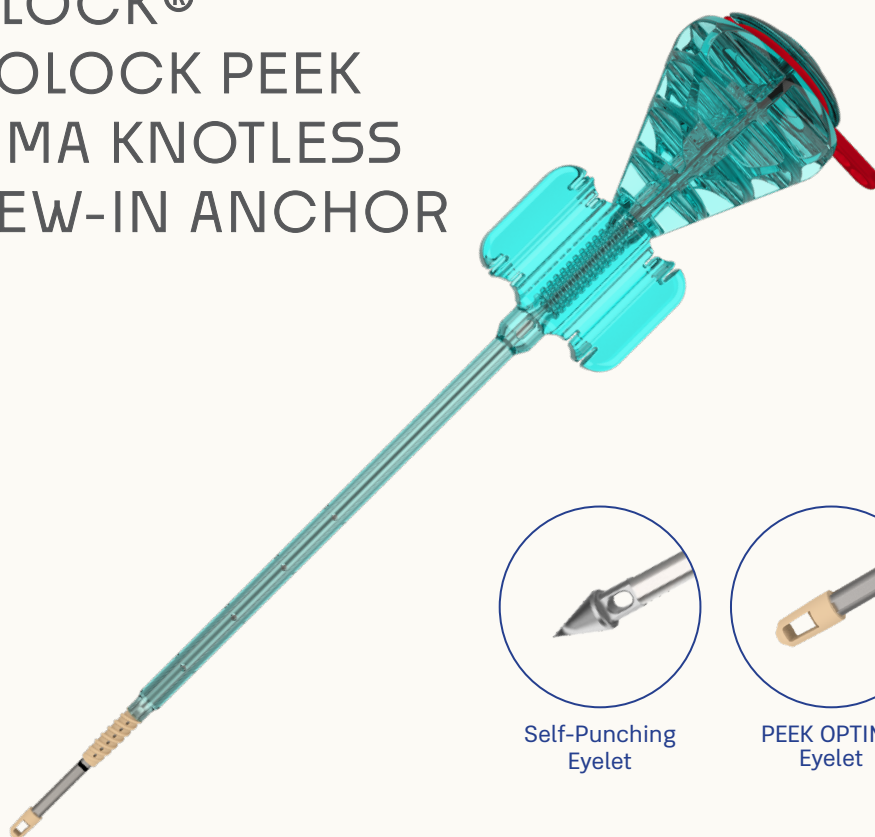

Self-Punching
Eyelet

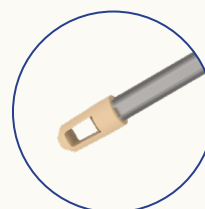

PEEK OPTIMA
Eyelet



PEEK OPTIMA
Anchor

GyroLock Knotless PEEK OPTIMA Screw-in Anchors consists of PEEK-OPTIMA (Poly-ether-ether-ketone) material. The GyroLock suture anchor comes with Dual start thread technology (fast inserting suture anchor) which allows anchor deployment in 3.5 rotation.

Vented and cannulated GyroLock anchors are designed for repair and reconstruction with soft-tissue graft. Vented design promotes bone growth. It is also available in a self-punching eyelet, designed for the lateral row in RC repair.

The GyroLock is designed for multiple procedures in single or double row for lateral row applications. GyroLock suture anchor eyelet can accommodate up to four 1.4mm tapes or four #2 sutures.

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-69	Awl/Tap for 4.75/5.5mm GyroLock Screw-In Anchor
7-013-70	Tap for 4.75/5.5mm GyroLock Screw-In Anchor

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SILICON COLOR
6-052-02	AUXILOCK® GyroLock 4.75mm x 22.5mm Knotless PEEK OPTIMA Screw-In Anchor	
6-052-03	AUXILOCK® Gyrolock 5.5mm x 22.5mm Knotless PEEK OPTIMA Screw-In Anchor	
6-052-04	AUXILOCK® GyroLock SP 4.75mm x 24.5mm Knotless PEEK OPTIMA Screw-In Anchor, Self-Punching	
6-052-05	AUXILOCK® GyroLock SP 5.5mm x 24.5mm Knotless PEEK OPTIMA Screw-In Anchor, Self-Punching	
6-052-11	AUXILOCK® Gyrolock 3.5mm x 15.8mm PEEK OPTIMA Screw-In Anchor With One 1.4mm Suture Tape: White/Blue, With Needle: MO-6	
6-052-12	AUXILOCK® Gyrolock 3.5mm x 15.8mm PEEK OPTIMA Screw-In Anchor With One 1.4mm Suture Tape: White/Black	
6-052-13	AUXILOCK® Gyrolock 3.5mm x 15.8mm PEEK OPTIMA Screw-In Anchor With One 2.0mm Suture Tape: Black/White	
6-052-14	AUXILOCK® Gyrolock 3.5mm x 15.8mm PEEK OPTIMA Screw-In Anchor With One 2.0mm Suture Tape: White/Blue, With Needle: MO-6	



AUXILOCK®
GYROLOCK PEEK
OPTIMA KNOTLESS
SCREW-IN ANCHOR



ORDERING INFORMATION



6-052-15	AUXILOCK® Gyrolock 4.75mm x 22.5mm PEEK OPTIMA Screw-In Anchor With One 2.0mm Suture Tape: White/Blue, With Needle: MO-6	 
6-052-16	AUXILOCK® Gyrolock 4.75mm x 22.5mm PEEK OPTIMA Screw-In Anchor With One 2.0mm Suture Tape: Black/White	 
6-052-17	AUXILOCK® Gyrolock 4.75mm x 22.5mm PEEK OPTIMA Screw-In Anchor With Two 2.0mm Suture Tape: White/Blue & Black/White, With Needle: MO-6	   
6-052-18	AUXILOCK® Gyrolock 4.75mm x 22.5mm PEEK OPTIMA Screw-In Anchor With Two 2.0mm Suture Tape: White/Blue & Black/White	   
6-052-19	AUXILOCK® Gyrolock 5.5mm x 22.5mm PEEK OPTIMA Screw-In Anchor With One 2.0mm Suture Tape: Black/White, With Needle: MO-6	 
6-052-20	AUXILOCK® Gyrolock 5.5mm x 22.5mm PEEK OPTIMA Screw-In Anchor With One 2.0mm Suture Tape: White/Blue	 
6-052-21	AUXILOCK® Gyrolock 5.5mm x 22.5mm PEEK OPTIMA Screw-In Anchor With Two 2.0mm Suture Tape: White/Blue & Black/White, With Needle: MO-6	   
6-052-22	AUXILOCK® Gyrolock 5.5mm x 22.5mm PEEK OPTIMA Screw-In Anchor With Two 2.0mm Suture Tape: White/Blue & Black/White	   

ORDERING INFORMATION


6-052-23	AUXILOCK® Gyrolock SP 4.75mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With One 2.0mm Suture Tape: White/Blue, With Needle: MO-6	 
6-052-24	AUXILOCK® Gyrolock SP 4.75mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With One 2.0mm Suture Tape: Black/White	 
6-052-25	AUXILOCK® Gyrolock SP 4.75mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With Two 2.0mm Suture Tape: White/Blue, Black/White With Needle: MO-6	   
6-052-26	AUXILOCK® Gyrolock SP 4.75mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With Two 2.0mm Suture Tape: White/Blue & Black/White	   
6-052-27	AUXILOCK® Gyrolock SP 5.5mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With One 2.0mm Suture Tape: Black/White, With Needle: MO-6	 
6-052-28	AUXILOCK® Gyrolock SP 5.5mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With One 2.0mm Suture Tape: White/Blue	 
6-052-29	AUXILOCK® Gyrolock SP 5.5mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With Two 2.0mm Suture Tape: White/Blue & Black/White, With Needle: MO-6	   
6-052-30	AUXILOCK® Gyrolock SP 5.5mm x 24.5mm PEEK OPTIMA Screw-In Anchor, Self-Punching With Two 2.0mm Suture Tape: White/Blue & Black/White	   



AUXILOCK® DUWLOC KNOTLESS PEEK OPTIMA PUSH-IN ANCHOR

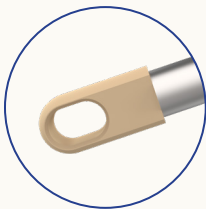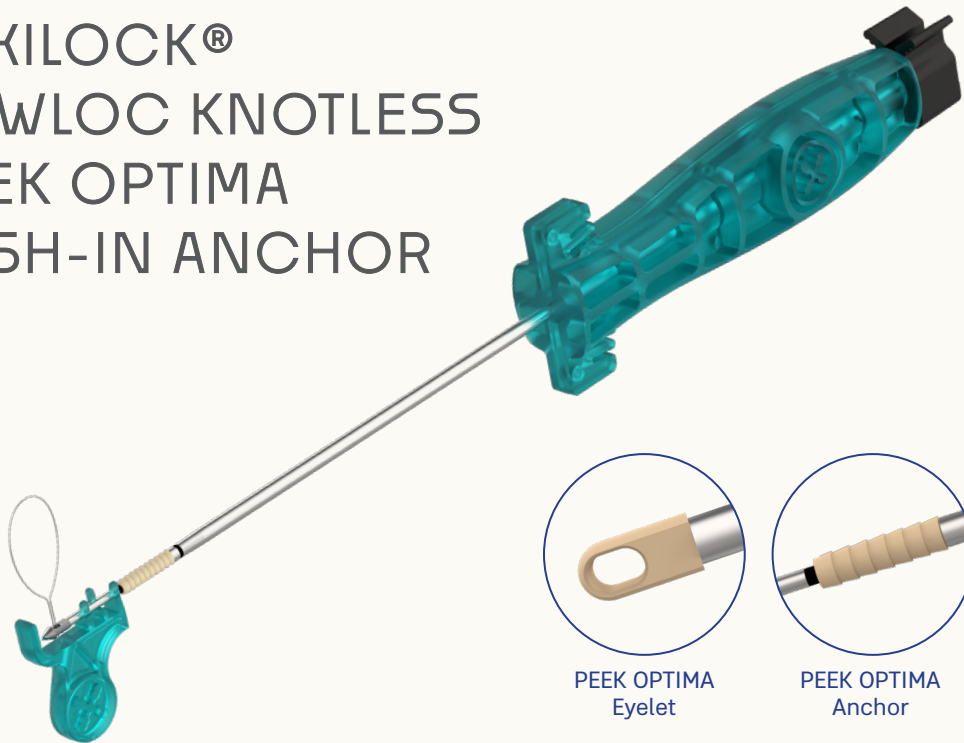

PEEK OPTIMA
Eyelet

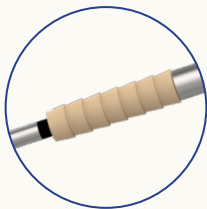

PEEK OPTIMA
Anchor



Self Punching
Eyelet

DUWLOC delivers secure, knotless rotator cuff repair with precise tension control-minimizing surgical steps and maximizing glenoid preservation.

AUXILOCK® DUWLOC PEEK OPTIMA Knotless Push-in Anchor is engineered for efficient, secure bankart and SLAP repair without the complications of arthroscopic knots. It is available in the diameter of 2.8mm & 3.5mm Its innovative “suture-first” technique allows surgeons to pass sutures through the labrum repair and fine-tune tissue tension before anchor insertion. It comes with a PEEK OPTIMA Eyelet a two piece design allowing for tension adjustment before locking. The DUWLOC anchor reduces the risk of knot-related impingement, offering multiple stitch configurations and precise tension control for optimal outcomes.

RECOMMENDED INSTRUMENTATION


CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-84	Ø3.2mm X 200mm Drill Bit	6-064-03
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-064-01, 6-064-03
7-013-02	8 Point Drill Guide, Length 100mm	6-064-01, 6-064-03
7-013-83	Ø2.5mm X 200mm Drill Bit	6-064-01

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-064-01	AUXILOCK® Duwloc 2.8mm x 15.5mm Knotless PEEK OPTIMA Push-In Anchor
6-064-02	AUXILOCK® Duwloc SP 2.8mm x 18.0mm Knotless PEEK OPTIMA Push-In Anchor
6-064-03	AUXILOCK® Duwloc 3.5mm x 20.0mm Knotless PEEK OPTIMA Push-In Anchor
6-064-04	AUXILOCK® Duwloc SP 3.5mm x 23.0mm Knotless PEEK OPTIMA Push-In Anchor
6-064-05	AUXILOCK® Duwloc 2.8mm x 15.5mm Knotless PEEK CF Push-In Anchor
6-064-06	AUXILOCK® Duwloc SP 2.8mm x 18.0mm Knotless PEEK CF Push-In Anchor
6-064-07	AUXILOCK® Duwloc 3.5mm x 20.0mm Knotless PEEK CF Push-In Anchor
6-064-08	AUXILOCK® Duwloc SP 3.5mm x 23.0mm Knotless PEEK CF Push-In Anchor
6-064-09	AUXILOCK® Duwloc 4.5mm x 27.0mm Knotless PEEK Optima Push-In Anchor
6-064-10	AUXILOCK® Duwloc SP 4.5mm x 30.0mm Knotless PEEK Optima Push-In Anchor
6-064-11	AUXILOCK® Duwloc 4.5mm x 27.0mm Knotless PEEK CF Push-In Anchor
6-064-12	AUXILOCK® Duwloc SP 4.5mm x 30.0mm Knotless PEEK CF Push-In Anchor
6-064-13	AUXILOCK® Duwloc 5.5mm x 27.0mm Knotless PEEK Optima Push-In Anchor
6-064-14	AUXILOCK® Duwloc SP 5.5mm x 31.5mm Knotless PEEK Optima Push-In Anchor
6-064-15	AUXILOCK® Duwloc 5.5mm x 27.0mm Knotless PEEK CF Push-In Anchor
6-064-16	AUXILOCK® Duwloc SP 5.5mm x 31.5mm Knotless PEEK CF Push-In Anchor



AUXILOCK® AUXPIK SUTURE- BASED ANCHOR

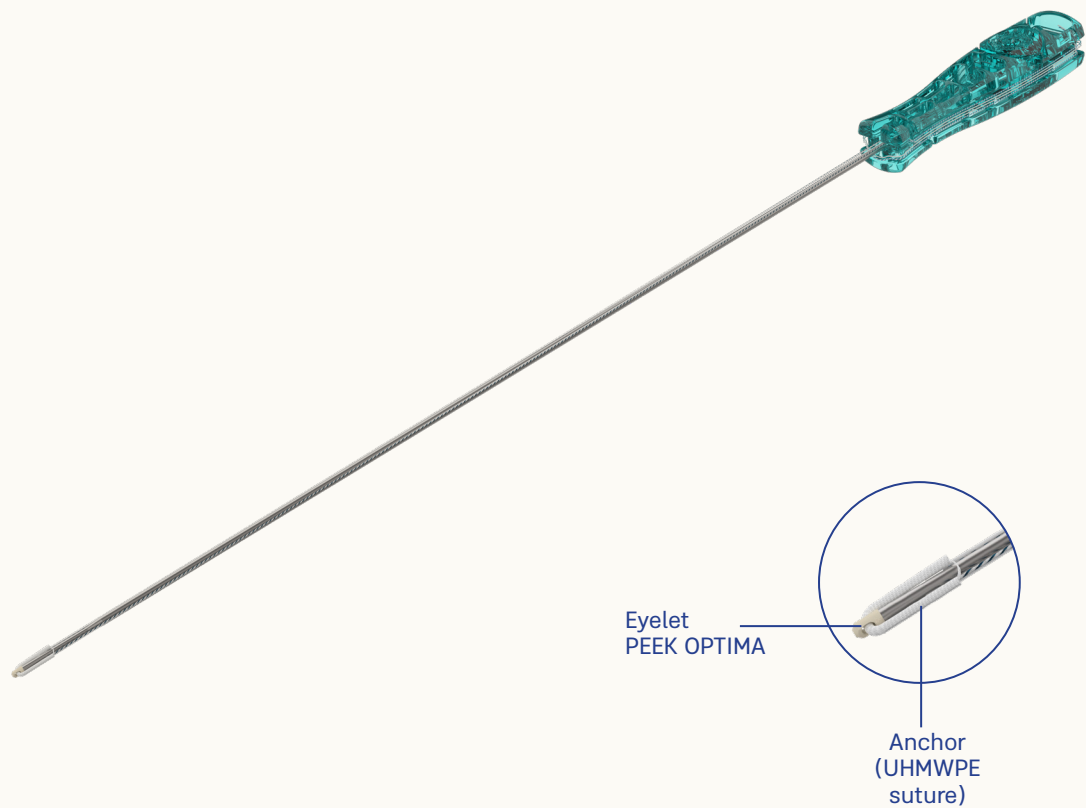

SHOULDER INSTABILITY REPAIR

AUXILOCK® AuxPik Suture-Based Anchor is made from UHMWPE suture anchor body with a PEEK OPTIMA eyelet tip. The AuxPik Suture-Based Anchor is designed for efficiency and ease of use. AuxPik Suture-Based Anchor is available in 1.8mm diameter with a BioBraid and suture tape. The anchor can be used for Bankart or SLAP repair surgeries. The anchor provides small foot print and subcortical fixation for anchor insertions. The drill guide and anchor driver combination are well designed to improve performance and reliability.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR
6-001-01	AUXILOCK® 1.8mm AuxPik Suture-Based Anchor With One #2 Bio-Braid: White/Blue	
6-001-08	AUXILOCK® 1.8mm AuxPik Suture-Based Anchor With One 1.4mm Suture Tape: White/Blue	

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-03	V-Mouth Drill Guide	6-001-01, 6-001-08
7-013-04	8 Point Drill Guide	6-001-01, 6-001-08
7-013-06	Trocar Tip Obturator	6-001-01, 6-001-08
7-013-53	Ø1.8mm X 325mm Drill Bit	6-001-01, 6-001-08



AUXILOCK® AUXSUTURE ANCHOR

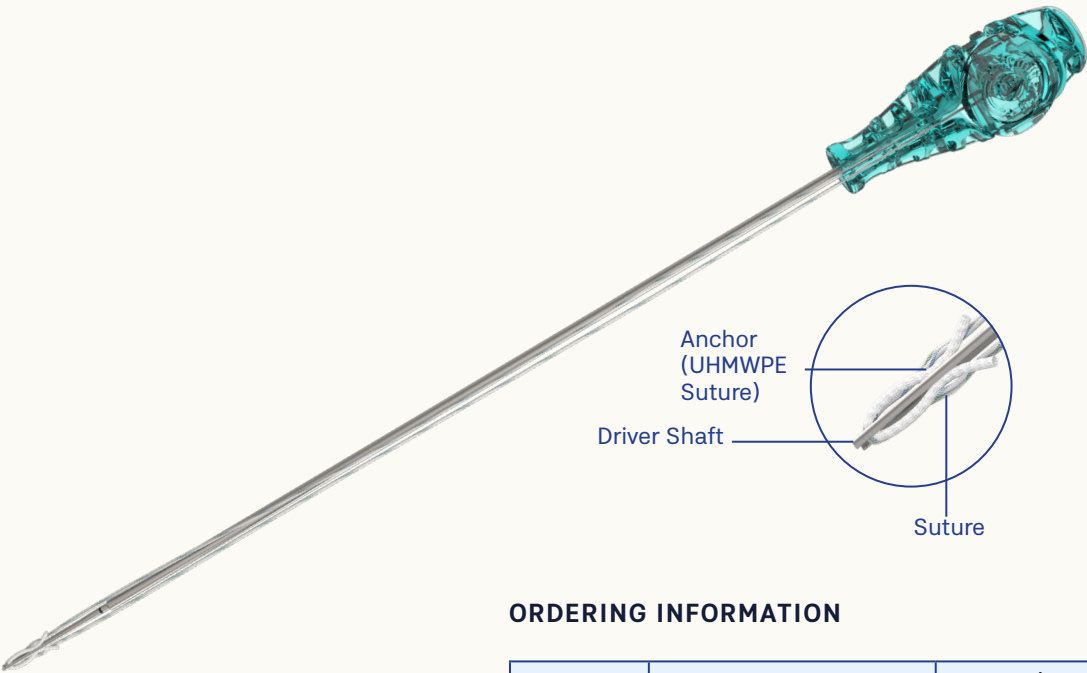


AUXILOCK® AuxSuture Anchor Provides numerous advantages over traditional anchors for shoulder instability repair procedures with its strong fixations, less bone removal and a push in technique. AuxSuture anchor is available in diameter of 1.5mm and 1.8mm.

It is available in single and double loaded options with various combination of BioBraid and suture tapes. The drill guide and 1.5mm drill bit are specially designed for AuxSuture Anchor while assuming an accurate hole placement for the anchor. The small 1.5mm footprint AuxSuture anchor provides placement and fixation advantages when bone real estate is limited, especially during revision surgeries. The anchor is also available with needles which are ideal for mini-open surgeries.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-034-01	AUXILOCK® 1.5mm AuxSuture Anchor With One #2 BioBraid: White/Blue		
6-034-02	AUXILOCK® 1.5mm AuxSuture Anchor With One 1.4mm Suture Tape: White/Blue		
6-034-03	AUXILOCK® 1.5mm Aux-Suture Anchor With One #2 BioBraid: White/Blue, With Needles: MO-6		
6-035-01	AUXILOCK® 1.8mm Aux-Suture Anchor With Two #2 BioBraid: White/Blue & White/Green		
6-035-02	AUXILOCK® 1.8mm AuxSuture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Black & White/Blue		
6-035-03	AUXILOCK® 1.8mm Aux-Suture Anchor With Two #2 BioBraid: White/Blue & White/Green, With Needles: MO-6		
6-035-04	AUXILOCK® 1.8mm AuxSuture Anchor With One #2 BioBraid and One 1.4mm Suture Tape: White/Black & White/Blue, With Needles: MO-6		

SHOULDER INSTABILITY REPAIR

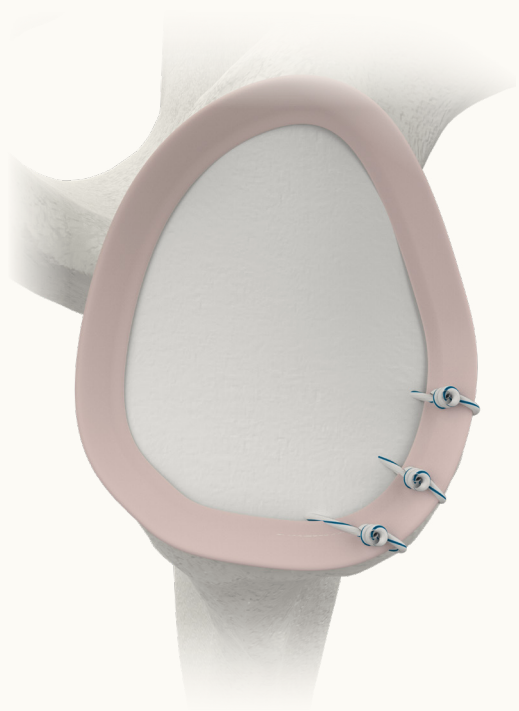


SHOULDER INSTABILITY REPAIR

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-03	V-Mouth Drill Guide	6-034-01, 6-034-02, 6-035-01, 6-035-02
7-013-04	8 Point Drill Guide	6-034-01, 6-034-02, 6-035-01, 6-035-02
7-013-05	V-Mouth Drill Guide - Curve	6-034-01, 6-034-02, 6-035-01, 6-035-02
7-013-06	Trocar Tip Obturator	6-034-01, 6-034-02, 6-035-01, 6-035-02
7-013-07	Trocar Tip Obturator - Curve	6-034-01, 6-034-02, 6-035-01, 6-035-02
7-013-20	Ø1.5mm, X 300mm Drill Bit, Stainless Steel	6-034-01, 6-034-02, 6-034-03, 6-035-01, 6-035-02, 6-035-03, 6-035-04
7-013-21	Ø1.5mm X 300mm Drill Bit, Nitinol	6-034-01, 6-034-02, 6-035-01, 6-035-02
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-034-03, 6-035-03, 6-035-04
7-013-02	8 Point Drill Guide, Length 100mm	6-034-03, 6-035-03, 6-035-04



AUXILOCK® PEEK CF PUSH-IN SUTURE ANCHOR

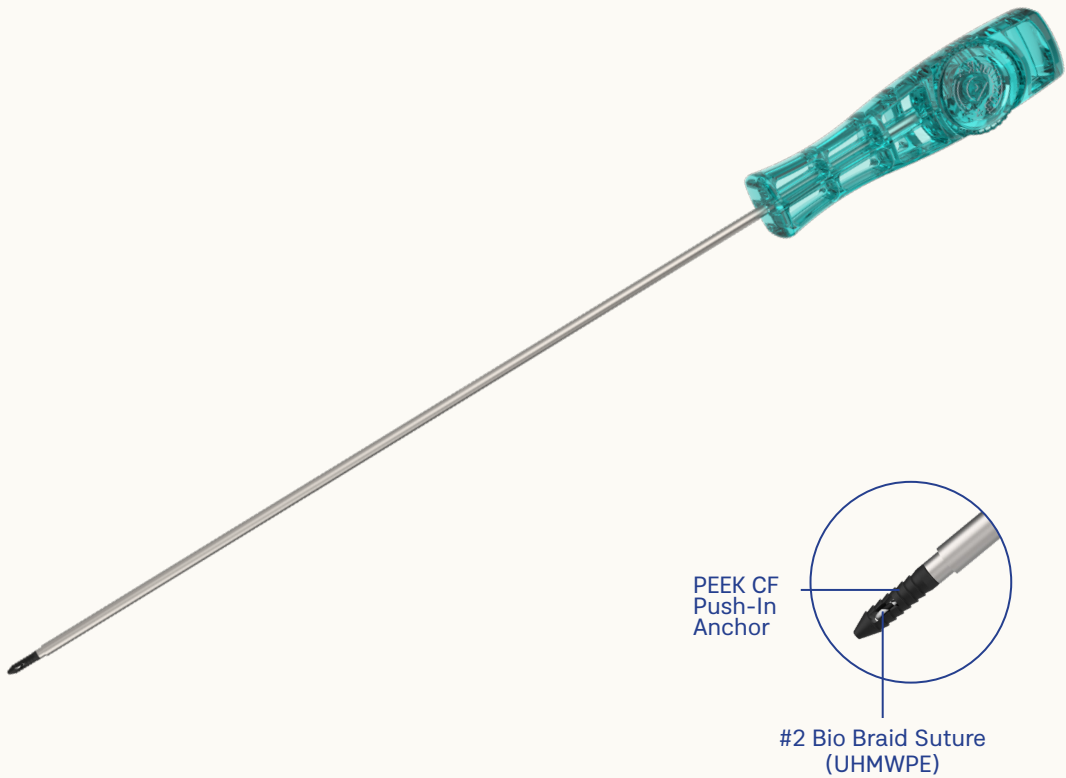

SHOULDER INSTABILITY REPAIR

AUXILOCK® PEEK CF Push-In Suture Anchor is a Push-in/ Bang-on anchor available in diameter of 2.8 and 3.5mm. It is available in single and double loaded sutures with and without needle options.

The PEEK CF Push-In Suture Anchors are made of carbon fiber-reinforced PEEK OPTIMA. The PEEK CF provides modulus of elasticity closely matching the cortical bone. It is non-absorbable, radiolucent, and MRI safe. The drill guide and drill bit are provided in the instrument set for the accurate placement of the anchor which minimizes the anchor slippage or breakage during the surgery. Anchor insertion and delivery are made simple by drilling a hole through a drill guide and inserting the anchor through the same drill guide into the drilled hole. The anchor is also available with needles which are ideal for mini-open Bankart or SLAP repair procedures.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-002-06	AUXILOCK® 2.8mm PEEK CF Push-In Suture Anchor With One #2 BioBraid: White/Blue		
6-002-07	AUXILOCK® 2.8mm PEEK CF Push-In Suture Anchor With One #2 Bio-Braid:White/Blue, With Needles: MO-6, Short Length		
6-002-01	AUXILOCK® 3.5mm PEEK CF Push-In Suture Anchor With Two #2 Bio-Braid:White/Blue & White/Black		
6-002-02	AUXILOCK® 3.5mm PEEK CF Push-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6,Short Length		




SHOULDER INSTABILITY REPAIR

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-03	V-Mouth Drill Guide	6-002-06, 6-002-01
7-013-04	8 Point Drill Guide	6-002-06, 6-002-01
7-013-06	Trocar Tip Obturator	6-002-06, 6-002-01
7-013-17	Ø2.4mm X 300mm Drill Bit	6-002-06
7-013-18	Ø2.8mm X 300mm Drill Bit	6-002-01
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-002-07, 6-002-02
7-013-12	Ø2.4mm X 200mm Drill Bit	6-002-07
7-013-11	Ø2.8mm X 200mm Drill Bit	6-002-07
7-013-02	8 Point Drill Guide, Length 100mm	6-002-07, 6-002-02



AUXILOCK® PEEK OPTIMA PUSH-IN SUTURE ANCHOR

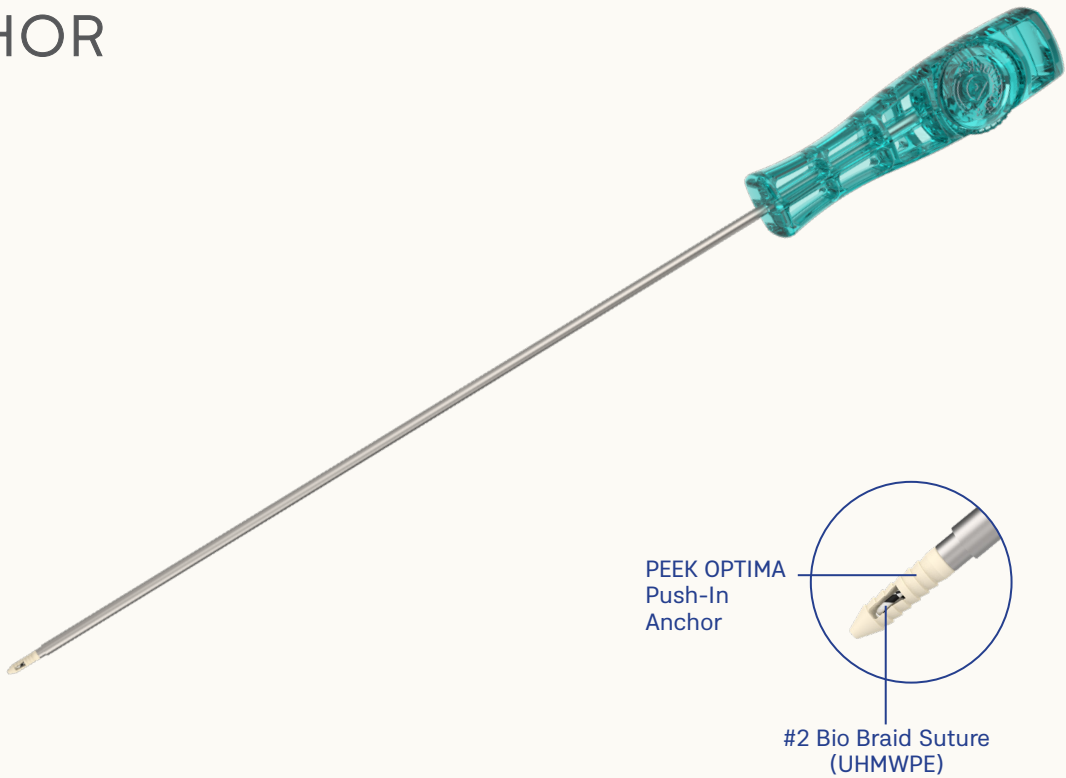

SHOULDER INSTABILITY REPAIR

AUXILOCK® PEEK OPTIMA Push-In Suture Anchor is a Push-in/ Bang-on anchor available in diameter of 2.8 and 3.5mm. It is available in single and double loaded sutures with and without needle options. The PEEK OPTIMA Push-In Suture Anchors are made of PEEK-OPTIMA (Poly-etheretherketone).

The PEEK OPTIMA Push-In Suture Anchors are non-absorbable PEEK suture anchors with a material eyelet. It provides superior abrasion resistance due to PEEK’s low coefficient of friction. It is considered to be Biocompatible, radiolucent, and MRI safe. The drill guide and the drill bit are provided in the instrument set for an accurate placement of the anchor which minimizes the anchor slippage or breakage during the surgery. Anchor insertion and delivery are made simple by drilling a hole through a drill guide and inserting the anchor through the same drill guide into the drilled hole.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-030-03	AUXILOCK 2.8mm PEEK OPTIMA Push-In Suture Anchor With One #2 BioBraid: White/Black		
6-030-04	AUXILOCK® 2.8mm PEEK OPTIMA Push-In Suture Anchor With One #2 Bio-Braid: White/Black, With Needles: MO-6, Short Length		
6-030-01	AUXILOCK® 3.5mm PEEK OPTIMA Push-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black		
6-030-02	AUXILOCK® 3.5mm PEEK OPTIMA Push-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6, Short Length		




SHOULDER INSTABILITY REPAIR

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-03	V-Mouth Drill Guide	6-030-03, 6-030-01
7-013-04	8 Point Drill Guide	6-030-03, 6-030-01
7-013-06	Trocar Tip Obturator	6-030-03, 6-030-01
7-013-17	Ø2.4mm X 300mm Drill Bit	6-030-03
7-013-18	Ø2.8mm X 300mm Drill Bit	6-030-01
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-030-04, 6-030-02
7-013-12	Ø2.4mm X 200mm Drill Bit	6-030-04
7-013-11	Ø2.8mm X 200mm Drill Bit	6-030-02
7-013-02	8 Point Drill Guide, Length 100mm	6-030-04



AUXILOCK® TITANIUM SCREW-IN SUTURE ANCHOR

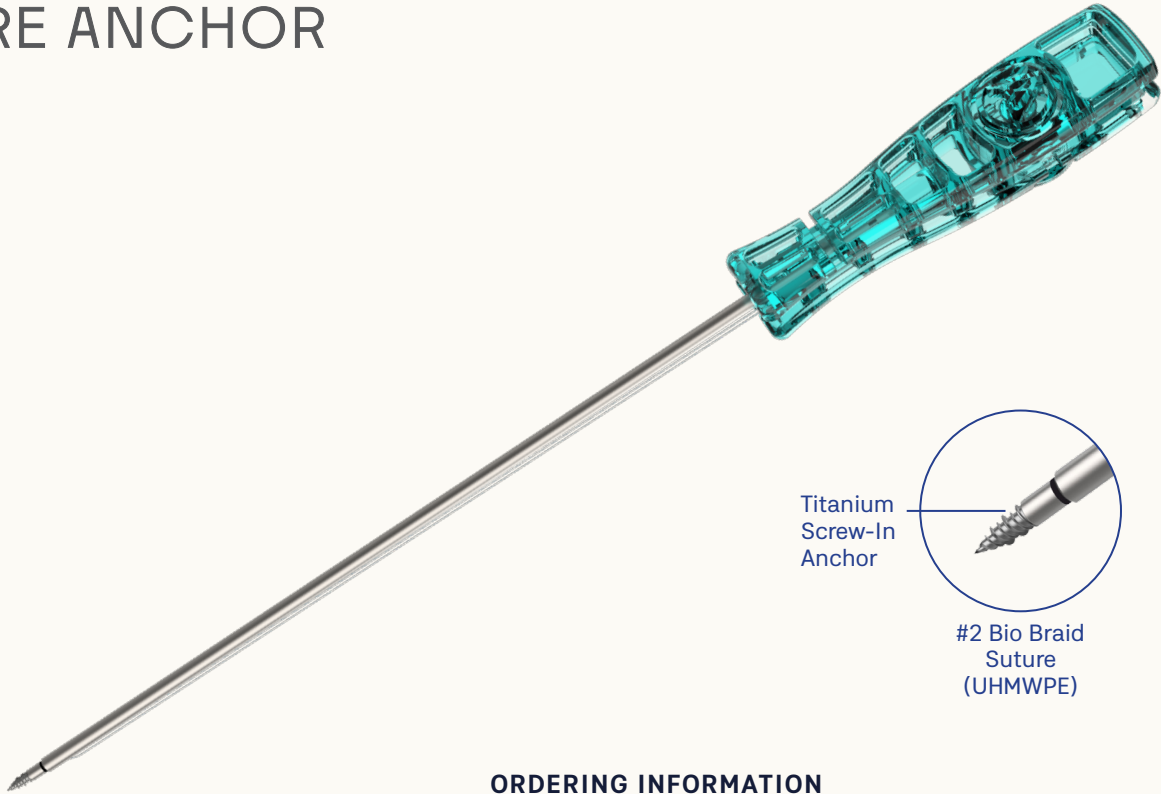


AUXILOCK® Titanium anchors are the screw-in anchors available in diameter of 2.8 and 3.5mm with single and double loaded options. It is a self-tapping suture anchor mostly used in Bankart, SLAP and elbow surgeries.

The anchor design incorporates a cortical thread to maximise pull-out strength in a cortical or hard bone. The anchor is designed for delivering its ultimate mechanical properties (pull-out strength, tensile strength, etc.) and ease-of-use. The drill guide and drill bit are provided in the instrument set for the accurate placement of the anchor. It in turn minimizes the anchor slippage or breakage during the surgery. The anchors are also available with needles, which are ideal for mini open elbow tendon repair and shoulder instability procedures.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION	SUTURE/TAPE COLOUR	NEEDLE
6-005-01	AUXILOCK® 2.8mm Titanium Screw-In Suture Anchor With One #2 BioBraid: White/Blue		
6-005-05	AUXILOCK® 2.8mm Titanium Screw-In Suture Anchor With One #2 BioBraid: White/Blue, With Needles: MO-6		
6-005-02	AUXILOCK® 3.5mm Titanium Screw-In Suture Anchor With Two #2 BioBraid: White/Blue & White/Black		
6-005-06	AUXILOCK® 3.5mm Titanium Screw-In Suture Anchor With Two #2 BioBraid: White & White/Blue, With Needles: MO-6		
6-005-13	AUXILOCK® 3.5mm Titanium Screw-In Suture Anchor With One #2 BioBraid: White/Blue, With Needles: MO-6		
6-005-14	AUXILOCK® 3.5mm Titanium Screw-In Suture Anchor With One #2 BioBraid: White/Blue		

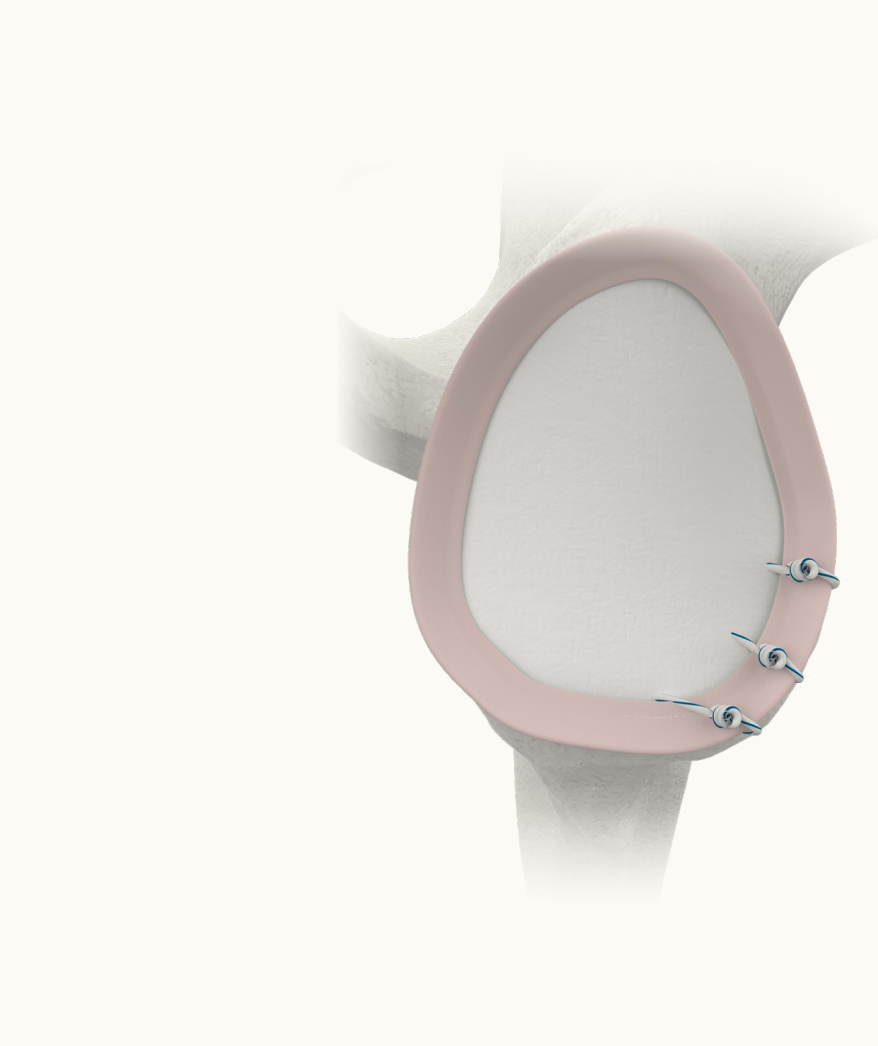


SHOULDER INSTABILITY REPAIR

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-03	V-Mouth Drill Guide	, 6-005-02, 6-005-14
7-013-04	8 Point Drill Guide	, 6-005-02, 6-005-14
7-013-06	Trocar Tip Obturator	, 6-005-02, 6-005-14
7-013-16	Ø1.8mm X 300mm Drill Bit	, 6-005-02, 6-005-14
7-013-02	8 Point Drill Guide, Length 100mm	6-005-05, 6-005-06, 6-005-13
7-013-10	Ø1.8mm X 200mm Drill Bit	6-005-05, 6-005-06, 6-005-13
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-005-05, 6-005-06, 6-005-13



AUXILOCK® KNOTLESS PEEK CF PUSH-IN SUTURE ANCHOR

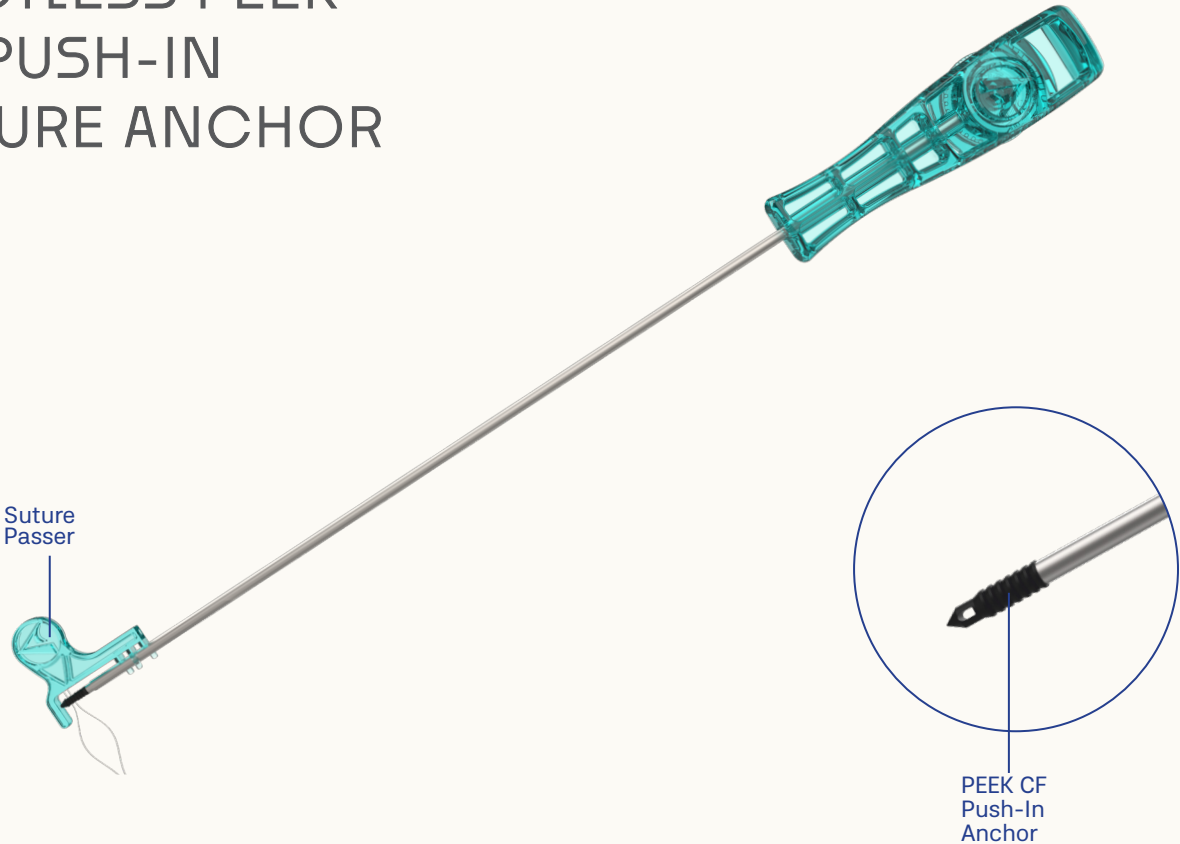

AUXILOCK® knotless push-in anchors are a step-saving alternative to conventional knotted suture anchors. AUXILOCK® knotless push-in anchors are recommended for use in Bankart and SLAP repair surgeries.

These Anchors are made of carbon fiber-reinforced PEEK-OPTIMA and are available in a variety of diameters of 2.8mm and 3.5mm providing intraoperative flexibility. The PEEK CF provides modulus of elasticity closely matching the cortical bone. It is non-absorbable, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION
6-003-01	AUXILOCK® 2.8mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer
6-003-02	AUXILOCK® 3.5mm Knotless PEEK CF Push-In Suture Anchor With Driver, With Suture Passer

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-003-01, 6-003-02
7-013-54	Ø2.5mm X 200 mm Drill Bit (For Knotless Anchor)	6-003-01
7-013-55	Ø3.2mm X 200 mm Drill Bit (For Knotless Anchor)	6-003-02
7-013-02	8 Point Drill Guide, Length 100mm	6-003-01, 6-003-02



AUXILOCK® KNOTLESS PEEK OPTIMA PUSH-IN SUTURE ANCHOR

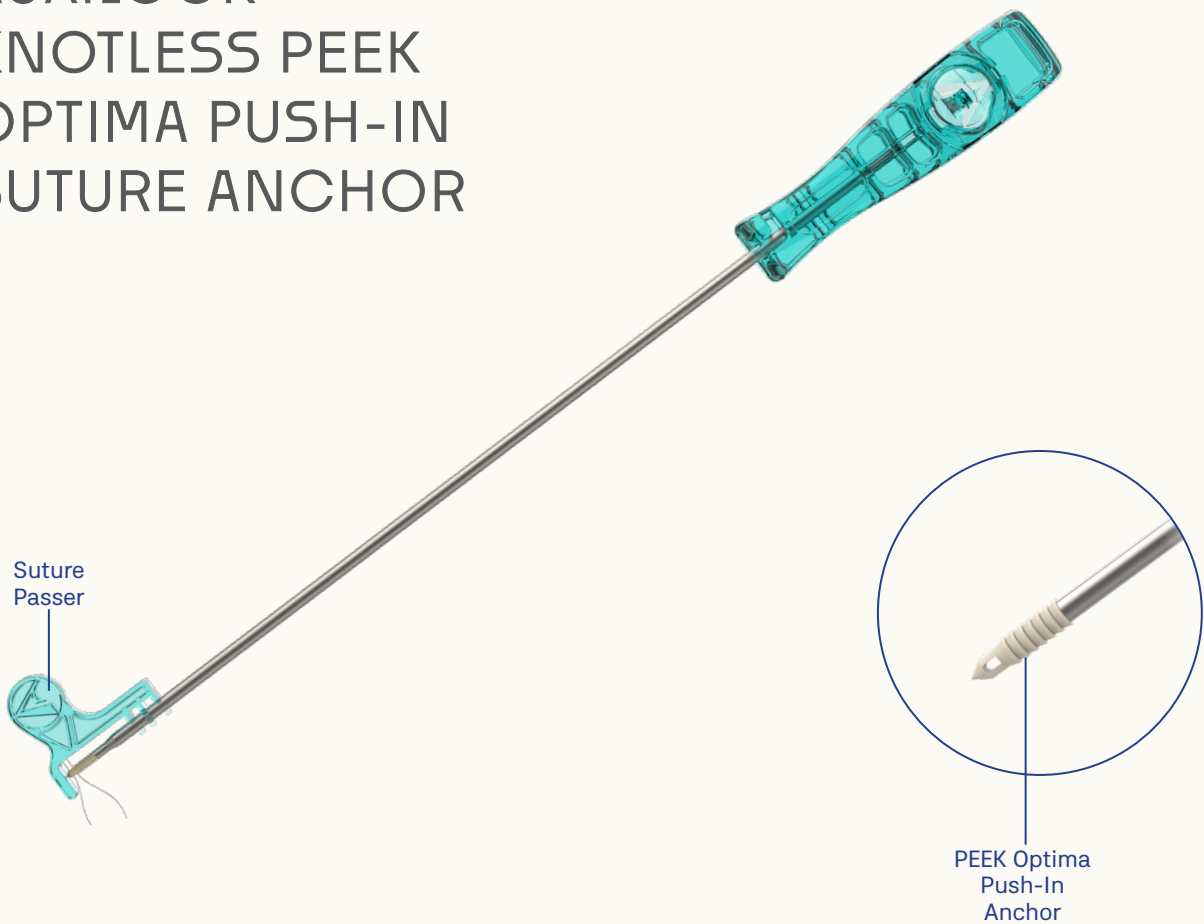

AUXILOCK® knotless push-in anchors are a step-saving alternative to conventional knotted suture anchors and are recommended for use in Bankart and SLAP repairs. The PEEK OPTIMA Push-In Suture Anchors are made of PEEK-OPTIMA (Poly-ether-ether-ketone).

These Anchors are available in a variety of diameters of 2.8mm and 3.5mm providing intraoperative flexibility. The PEEK OPTIMA Push-In Suture Anchor is a one piece body anchor with eyelet. It provides superior abrasion resistance due to PEEK's low coefficient friction. It is Biocompatible, radiolucent, and MRI safe. The knotless technology of Push-in anchor also eliminates knot stacks associated with soft tissue irritation.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-031-01	AUXILOCK® 2.8mm Knotless PEEK OPTIMA Push-In Suture Anchor With Driver, With Suture Passer
6-031-02	AUXILOCK® 3.5mm Knotless PEEK OPTIMA Push-In Suture Anchor With Driver, With Suture Passer

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION	COMPATIBLE IMPLANTS
7-013-01	8 Point Slotted Drill Guide, Length 190mm	6-031-01, 6-031-02
7-013-54	Ø2.5mm X 200 mm Drill Bit (For Knotless Anchor)	6-031-01
7-013-55	Ø3.2mm X 200 mm Drill Bit (For Knotless Anchor)	6-031-02
7-013-02	8 Point Drill Guide, Length 100mm	6-031-01, 6-031-02



AUXILOCK® GYROLOCK PEEK OPTIMA KNOTLESS SCREW-IN ANCHOR



GyroLock Knotless PEEK OPTIMA Screw-in Anchors consists of PEEK-OPTIMA (Poly-ether-ether-ke-tone) material. The GyroLock suture anchor comes with Dual start thread technology (fast inserting suture anchor) which allows anchor deployment in 3.5 rotation.

Vented and cannulated GyroLock anchors are designed for repair and reconstruction with soft-tissue graft. GyroLock suture anchor eyelet can accommodate up to two 1.4mm tapes or two #2 sutures. Vented design promotes bone growth.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION
6-052-01	AUXILOCK® GyroLock 3.5mm x 15.8mm Knotless PEEK OPTIMA Screw-In Anchor

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-71	Tap for 3.5mm GyroLock Screw-In Anchor



AUXILOCK® TENDO PEEK CF SCREW-IN ANCHOR

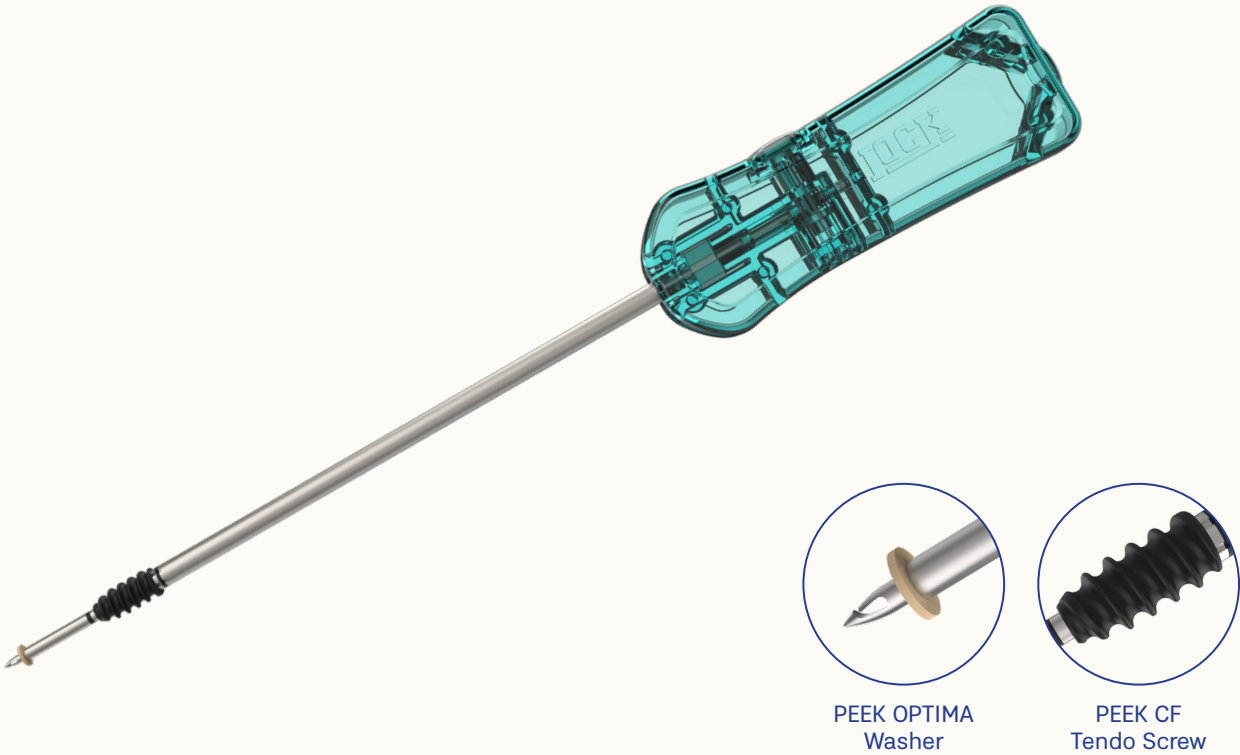

AUXILOCK® Tendo PEEK CF Screw-In Anchor is ideal for proximal and distal biceps tenodesis repair. It is available in three diameters 7.0, 8.0 & 9.0mm & two lengths 10mm and 15mm, with PEEK CF interference/tenodesis screws and a PEEK OPTIMA washer.

The Eyelet is specifically designed to hold the biceps tendon at the bottom of the drill hole. After fixing up the tendon into the drill hole, it can be locked into the bone by advancing the tenodesis screw. The anchor is designed to save steps and minimize the length of the procedure.

ORDERING INFORMATION


CODE	PRODUCT DESCRIPTION
6-043-01	AUXILOCK® Tendo 7.0mm x 15mm PEEK CF Screw-In Anchor
6-043-02	AUXILOCK® Tendo 8.0mm x 15mm PEEK CF Screw-In Anchor
6-043-03	AUXILOCK® Tendo 9.0mm x 15mm PEEK CF Screw-In Anchor
6-043-04	AUXILOCK® Tendo 7.0mm x 10mm PEEK CF Screw-In Anchor
6-043-05	AUXILOCK® Tendo 8.0mm x 10mm PEEK CF Screw-In Anchor
6-043-06	AUXILOCK® Tendo 9.0mm x 10mm PEEK CF Screw-In Anchor

RECOMMENDED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-013-58	Graft Sizer, 6.0mm to 9.0mm (1.0mm Steps)
7-013-59	Ø6.0mm Trocar Tip Headed Reamer
7-013-60	Ø7.0mm Trocar Tip Headed Reamer
7-013-61	Ø8.0mm Trocar Tip Headed Reamer
7-013-62	Ø9.0mm Trocar Tip Headed Reamer



AUXILOCK® TENDO PEEK OPTIMA SCREW-IN ANCHOR

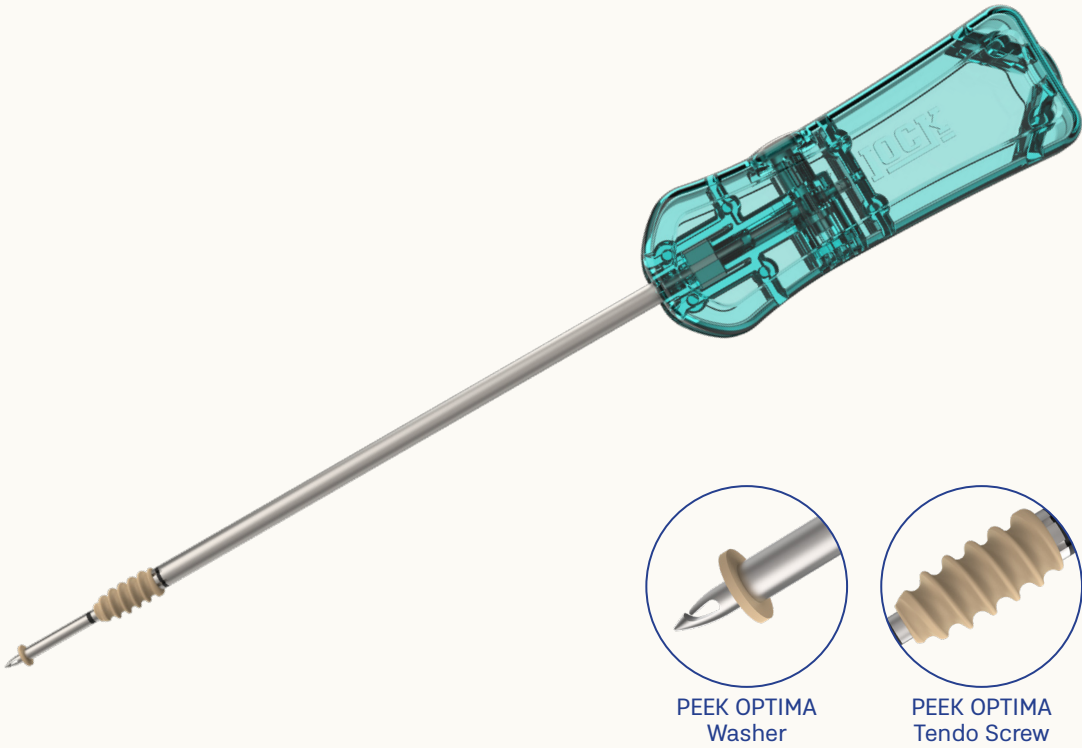

BICEPS TENODESIS REPAIR

AUXILOCK® Tendo PEEK OPTIMA Screw-In Anchor is ideal for proximal and distal biceps tenodesis repair. It is available in three diameters 7.0, 8.0 & 9.0mm & two lengths 10mm and 15mm, with PEEK OPTIMA interference/tenodesis screws and a PEEK OPTIMA washer.

The Eyelet is specifically designed to hold the biceps tendon at the bottom of the drill hole. After fixing up the tendon into the drill hole, it can be locked into the bone by advancing the tenodesis screw. The anchor is designed to save steps and minimize the length of the procedure.

ORDERING INFORMATION



CODE	PRODUCT DESCRIPTION
6-042-01	AUXILOCK® Tendo 7.0mm x 15mm PEEK OPTIMA Screw-In Anchor
6-042-02	AUXILOCK® Tendo 8.0mm x 15mm PEEK OPTIMA Screw-In Anchor
6-042-03	AUXILOCK® Tendo 9.0mm x 15mm PEEK OPTIMA Screw-In Anchor
6-042-04	AUXILOCK® Tendo 7.0mm x 10mm PEEK OPTIMA Screw-In Anchor
6-042-05	AUXILOCK® Tendo 8.0mm x 10mm PEEK OPTIMA Screw-In Anchor
6-042-06	AUXILOCK® Tendo 9.0mm x 10mm PEEK OPTIMA Screw-In Anchor

RECOMMENDED INSTRUMENTATION


CODE	PRODUCT DESCRIPTION
7-013-58	Graft Sizer, 6.0mm to 9.0mm (1.0mm Steps)
7-013-59	Ø6.0mm Trocar Tip Headed Reamer
7-013-60	Ø7.0mm Trocar Tip Headed Reamer
7-013-61	Ø8.0mm Trocar Tip Headed Reamer
7-013-62	Ø9.0mm Trocar Tip Headed Reamer



AUXILOCK® AC JOINT KIT

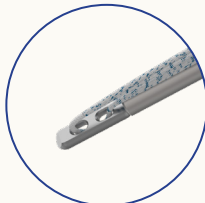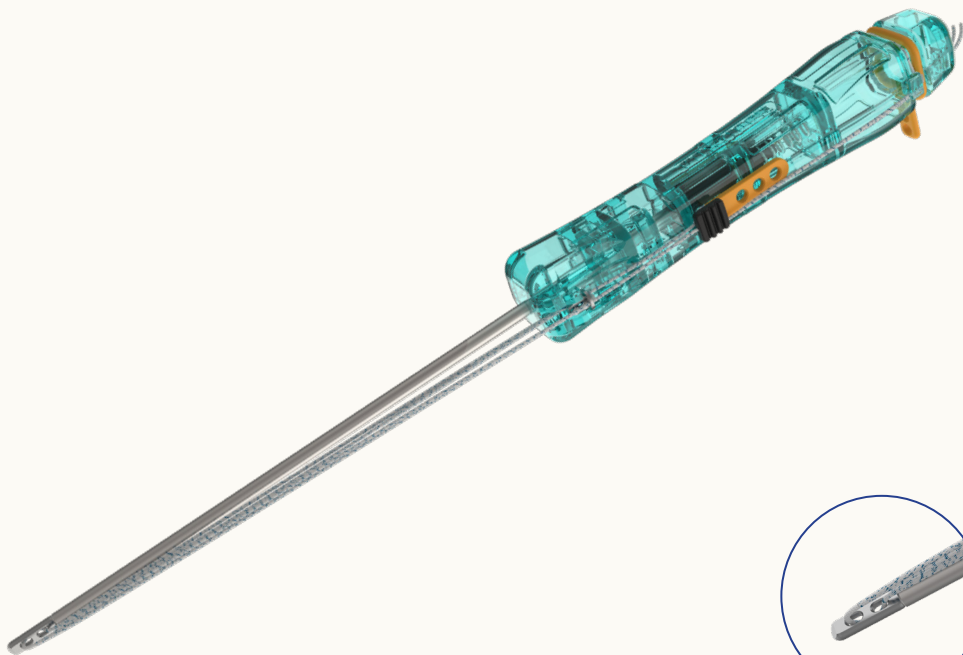

Oblong Shaped
Button (L-10mm
X W-3.2mm X
T-1.5mm).

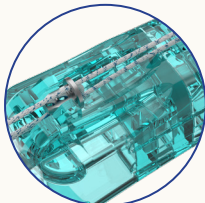

Round Button
5.5mm

AUXILOCK® AC Joint Kit can be used for Acromioclavicular Joint (AC Joint) separation repair. The kit comprises of adjustable loop with round shaped button and oblong shaped button. A tunnel is made through the clavicle to coracoid using 1.5mm K-Wire and 3.5mm Cannulated Drill Bit.

The kit is passed through a 3.5mm drilling tunnel from clavicle to coracoid, providing fixation with a round button on the clavicle and an oblong-shaped button at the inferior side of the coracoid. It is ideal for Rockwood type III, IV & V classification. The AC Joint Kit provides a double locking mechanism which eliminates the need for knot tying and also provides very strong compression on the clavicles.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-041-02	AUXILOCK® AC Joint Kit


RELATED INSTRUMENTATION (INCLUDE IN KIT)

CODE	PRODUCT DESCRIPTION
7-147-01	Cannulated Drill Bit with Quick Coupling End, Ø3.5mm x Length 175mm
7-147-02	Ø1.5mm K-Wire



AUXILOCK® AUXFIX BUTTON

Slotted Holes
to load Suture
Tapes

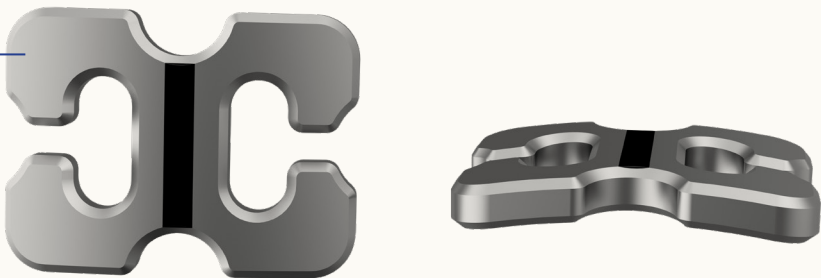


AUXILOCK® AuxFix Button is a titanium button with holes that allow use of multiple suture tapes for AC joint reduction. It provides very strong compression on the clavicles. The button is laid on the outer surface of the clavicle and underneath the coracoid. Only the suture material is passed through the clavicle and coracoid tunnel ensuring minimal bone loss during the procedures. Designed for use in chronic AC joint cases.

ORDERING INFORMATION

CODE	PRODUCT DESCRIPTION
6-037-01	AUXILOCK® AuxFix Button (AC Joint)

RELATED INSTRUMENTATION

CODE	PRODUCT DESCRIPTION
7-014-68	Passing Pin Ø2.4mm X 250mm
7-013-57	Ø4.5mm Cannulated Headed Reamer

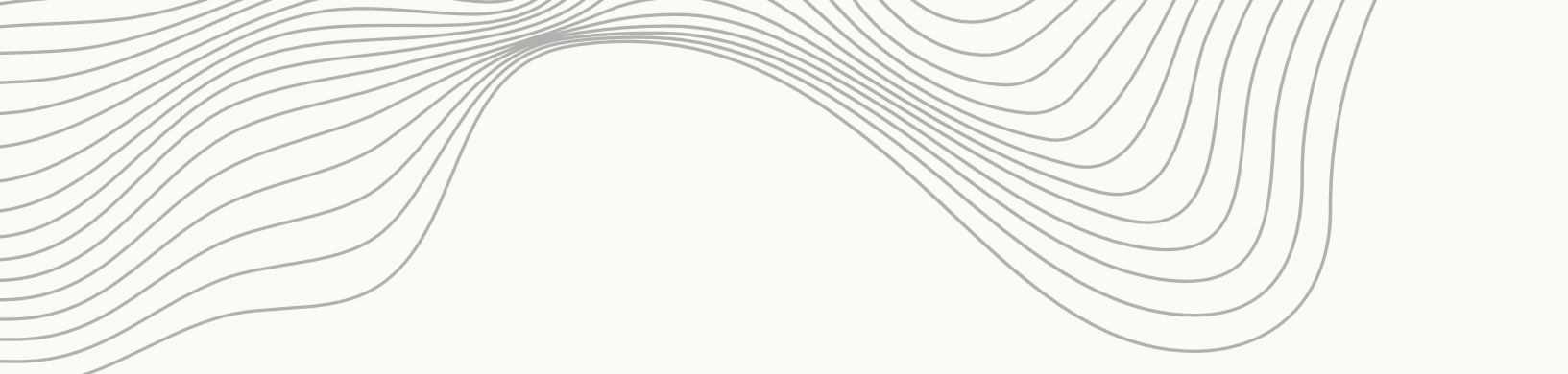

INSTRUMENTS



7-013-01

8 Point Slotted Drill Guide, Length 190mm



7-013-02

8 Point Drill Guide, Length 100mm



7-013-03

V-Mouth Drill Guide



7-013-04

8 Point Drill Guide



7-013-05

V-Mouth Drill Guide - Curve



7-013-06

Trocar Tip Obturator






7-013-07

Trocac Tip Obturator - Curve




7-013-10

Ø1.8mm X 200mm Drill Bit (For 2.8mm & 3.5mm Ti Needle Anchors)

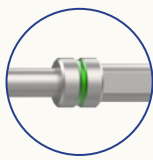

7-013-11

Ø2.8mm X 200mm Drill Bit (For 3.5mm Push-In Needle Anchor)

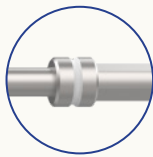

7-013-12

Ø2.4mm X 200mm Drill Bit (For 2.8mm Push-In Needle Anchor)




7-013-13

Ø3.2mm X 200mm Drill Bit (For 3.2mm AuxPix Anchor)



7-013-16


Ø1.8mm X 300mm Drill Bit (For 2.8mm & 3.5mm Ti Anchors)

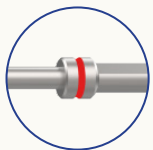




7-013-17

Ø2.4mm X 300mm Drill Bit (For 2.8mm Push-In Anchor)

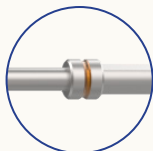


7-013-18

Ø2.8mm X 300mm Drill Bit (For 3.5mm Push-In Anchor)

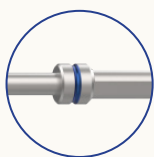

7-013-20

Ø1.5mm X 300mm Drill Bit, Stainless Steel (For 1.5mm & 1.8mm AuxSuture Anchor)

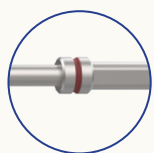


7-013-21

Ø1.5mm X 300mm Drill Bit, Nitinol (For 1.5mm & 1.8mm AuxSuture Anchor - Curve)




7-013-53

Ø1.8mm X 325 mm Drill Bit (For 1.8mm AuxPix Anchor)

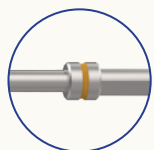

7-013-54

Ø2.5mm X 200 mm Drill Bit (For 2.8mm Knotless Anchor)

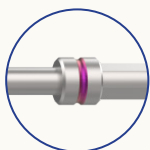




7-013-55

Ø3.2mm X 200 mm Drill Bit (For 3.5mm Knotless Anchor)



7-013-33

Switching Stick

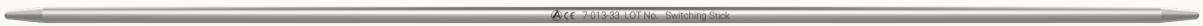

7-013-22

Ø2.9mm Awl (For 2.9mm AuxSuture Anchor)



7-013-23

Ø3.2mm Awl (For 3.2mm AuxPix Anchor & 4.5mm Knotless Anchor)



7-013-24

Ø4.1mm Awl (For Rotador & 5.5mm Knotless Anchor)



7-013-25

Awl For 4.5/5.5mm PEEK Screw-In Anchor (For PEEK Screw-In Anchor)





7-013-26

Awl For 6.5mm PEEK Screw-In Anchor (For PEEK Screw-In Anchor)



7-013-27

Awl/ Tap For 4.5/5.5mm PEEK Screw-In Anchor (For PEEK Screw-In Anchor)



7-013-28

Awl/ Tap For 6.5mm PEEK Screw-In Anchor (For PEEK Screw-In Anchor)



7-013-29

Tap For 4.5/5.5mm PEEK Screw-In Anchor (For PEEK Screw-In Anchor)



7-013-30

Tap For 6.5mm PEEK Screw-In Anchor (For PEEK Screw-In Anchor)



7-013-32

4.5/5.5mm Knotles PEEK Anchor Driver





7-013-34

Arthroscopic Probe - Shoulder



7-013-35

Tissue Elevator, 15° Up



7-013-36

Tissue Elevator, 15° Down



7-013-37

Rasp, 15° Up




7-013-38

Rasp, 15° Down



7-013-39

Ring Curette

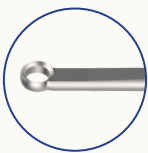




7-013-40

Knot Pusher



7-013-44

Hammer

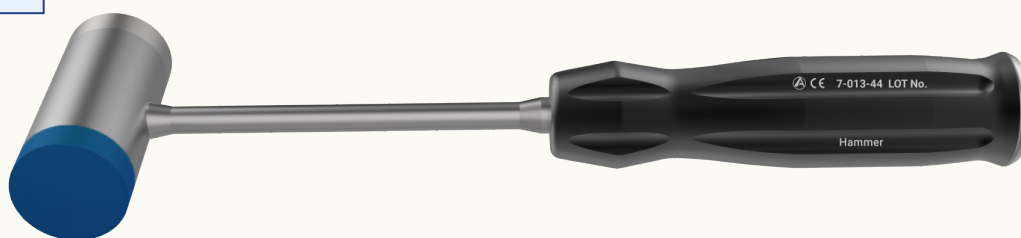


7-013-57

Ø4.5mm Cannulated Headed Reamer



7-013-58

Graft Sizer, 6.0mm to 9.0mm (1.0mm Steps)

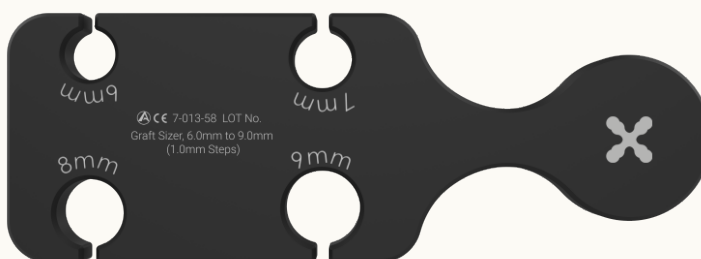

CODE	DIAMETER
7-013-59	Ø6.0mm Trocar Tip Headed Reamer
7-013-60	Ø7.0mm Trocar Tip Headed Reamer
7-013-61	Ø8.0mm Trocar Tip Headed Reamer
7-013-62	Ø9.0mm Trocar Tip Headed Reamer

Tocar Tip Headed Reamer



7-014-68

Passing Pin, Ø2.4mm x 250mm





7-013-69

Awl/Tap for 4.75/5.5mm Gyrolock Screw-In Anchor



7-013-70

Tap for 4.75/5.5mm Gyrolock Screw-In Anchor



7-013-71

Tap for 3.5mm Gyrolock Screw-In Anchor



7-013-83

Ø2.5mm X 200mm Drill Bit for 2.8mm (Dowloc Anchor)



7-013-84

Ø3.2mm X 200mm Drill Bit for 3.5mm (Dowloc Anchor)





OPTIONAL INSTRUMENTATION

7-013-63

Tissue Elevator, 30° Up



7-013-64

Tissue Elevator, 30° Down



Reusable Suture Passer



REUSABLE SUTURE PASSER

CODE	PRODUCT DESCRIPTION	HANDLE COLOR
7-013-75	Reusable Suture Passer, Straight	Black
7-013-76	Reusable Suture Passer, 25° Left	Red
7-013-77	Reusable Suture Passer, 45° Left	Purple
7-013-78	Reusable Suture Passer, 25° Right	Teal
7-013-79	Reusable Suture Passer, 45° Right	Green



Straight




25° Left



25° Right



45° Left

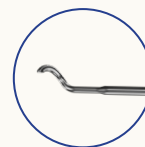


45° Right



7-013-80

AVAILABLE IN STERILE EO


Nitinol Loop for Reusable Suture Passer

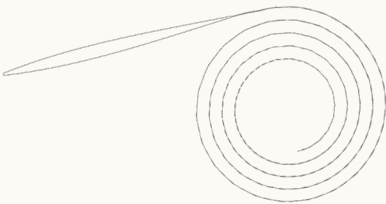

7-013-81

AVAILABLE IN STERILE EO

FlexPass Lasso

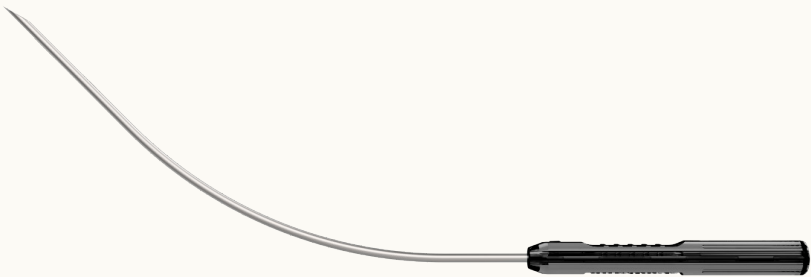

7-013-82

Blunt Tip Obturator



OPTIONAL INSTRUMENTATION

CODE	SET CONSISTING OF	UNIT
7-013-63	Tissue Elevator, 30° Up	1
7-013-64	Tissue Elevator, 30° Down	1
7-013-75	Reusable Suture Passer, Straight	1
7-013-76	Reusable Suture Passer, 25 °Left	1
7-013-77	Reusable Suture Passer, 45 °Left	1
7-013-78	Reusable Suture Passer, 25 °Right	1
7-013-79	Reusable Suture Passer, 45 °Right	1
7-013-80	Nitinol Loop for Reusable Suture Passer	1
7-013-81	FlexPass Lasso	1
7-013-82	Blunt Tip Obturator	1




General Hand Instruments of Shoulder

7-013-41




Suture Retriever

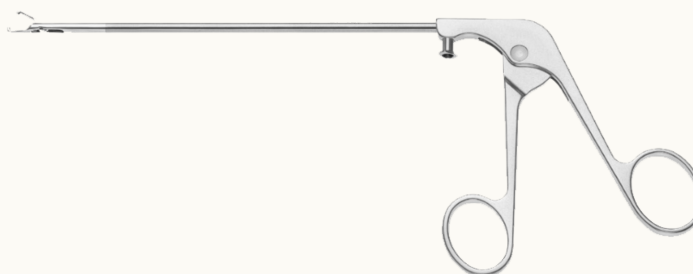

7-013-42




Tissue Grasper

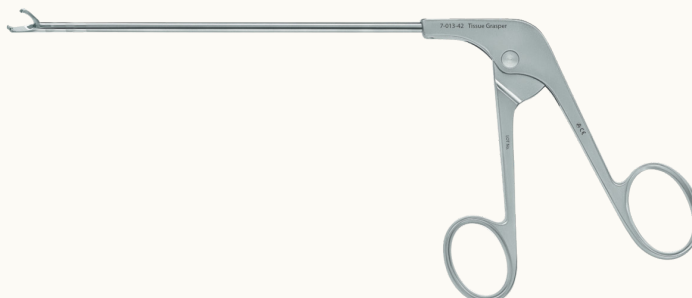

7-013-43




Suture Cutter/Knot Pusher With Lock, Ø4.2mm,
V-Knotch, Left Opened

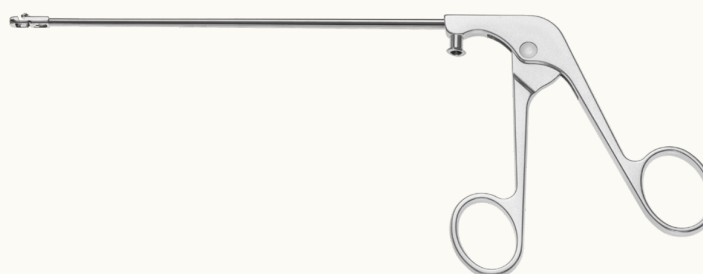

7-013-47



ESCORPIAN Suture Passer

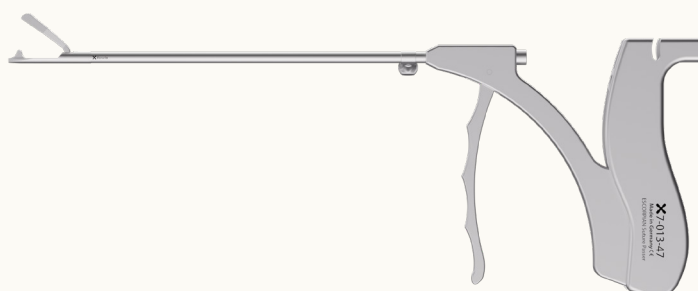



7-013-48

Needle With Lateral Opening For ESCORPIAN Suture Passer

AVAILABLE IN STERILE EO



CODE	DIAMETER
7-013-49	Smart Hook Suture passer, Straight
7-013-50	Smart Hook Suture passer, Right
7-013-51	Smart Hook Suture passer, Left

AVAILABLE IN STERILE EO

Smart Hook Suture passer



7-182-01


Container for Shoulder Arthroscopy General instruments





7-182-02

Lid for Shoulder Arthroscopy General Instruments

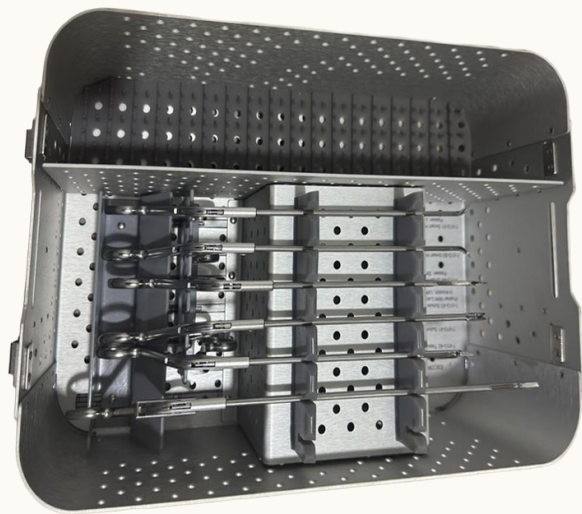

General Hand Instruments of Shoulder

CODE	SET CONSISTING OF	UNIT
7-013-41	Suture Retriever	1
7-013-42	Tissue Grasper	1
7-013-43	Suture Cutter/Knot Pusher With Lock, Ø4.2mm, V-Knotch, Left Opened	1
7-013-47	ESCORPIAN Suture Passer	1
7-013-48	Needle With Lateral Opening For ESCORPIAN Suture Passer	1
7-013-49	Smart Hook Suture passer, Straight	1
7-013-50	Smart Hook Suture passer, Right	1
7-013-51	Smart Hook Suture passer, Left	1
7-182-01	Container for Shoulder Arthroscopy General instruments	1
7-182-02	Lid for Shoulder Arthroscopy General instruments	1





ORDER INFORMATION

7-013	SHOULDER ARTHROSCOPY INSTRUMENT SET
-------	-------------------------------------

CODE	SET CONSISTING OF	UNIT
7-013-01	8 Point Slotted Drill Guide, Length 190mm	1
7-013-02	8 Point Drill Guide, Length 100mm	1
7-013-03	V-Mouth Drill Guide	1
7-013-04	8 Point Drill Guide	1
7-013-05	V-Mouth Drill Guide - Curve	1
7-013-06	Trocar Tip Obturator	1
7-013-07	Trocar Tip Obturator - Curve	1
7-013-10	Ø1.8mm X 200mm Drill Bit	1
7-013-11	Ø2.8mm X 200mm Drill Bit	1
7-013-12	Ø2.4mm X 200mm rill Bit(For 2.8mm Push-In Needle Anchor)	1
7-013-13	Ø3.2mm X 200mm Drill Bit (For 3.2mm AuxPix Anchor)	1
7-013-16	Ø1.8mm X 300mm Drill Bit	1
7-013-17	Ø2.4mm X 300mm Drill Bit	1
7-013-18	Ø2.8mm X 300mm Drill Bit	1
7-013-20	Ø1.5mm X 300mm Drill Bit, Stainless Steel	1
7-013-21	Ø1.5mm X 300mm Drill Bit, Nitinol	1
7-013-22	Ø2.9mm Awl	1
7-013-23	Ø3.2mm Awl	1
7-013-24	Ø4.1mm Awl	1
7-013-25	Awl For 4.5/5.5mm PEEK Screw-In Anchor	1
7-013-26	Awl For 6.5mm PEEK Screw-In Anchor	1
7-013-27	Awl/ Tap For 4.5/5.5mm PEEK Screw-In Anchor	1
7-013-28	Awl/ Tap For 6.5mm PEEK Screw-In Anchor	1
7-013-29	Tap For 4.5/5.5mm PEEK Screw-In Anchor	1
7-013-30	Tap For 6.5mm PEEK Screw-In Anchor	1
7-013-32	4.5/5.5mm Knotless PEEK Anchor Driver	1
7-013-33	Switching Stick	1
7-013-34	Arthroscopic Probe - Shoulder	1
7-013-35	Tissue Elevator, 15° Up	1
7-013-36	Tissue Elevator, 15° Down	1
7-013-37	Rasp, 15° Up	1
7-013-38	Rasp, 15° Down	1
7-013-39	Ring Curette	1
7-013-40	Knot Pusher	1
7-013-44	Hammer	1
7-013-53	Ø1.8mm X 325 mm Drill Bit	1
7-013-54	Ø2.5mm X 200 mm Drill Bit (For Knotless Anchor)	1
7-013-55	Ø3.2mm X 200 mm Drill Bit (For Knotless Anchor)	1
7-014-68	Passing Pin, Ø2.4mm x 250mm	1
7-013-57	Ø4.5mm Cannulated Headed Reamer	1



CODE	SET CONSISTING OF	UNIT
7-013-58	Graft Sizer, 6.0mm to 9.0mm (1.0mm Steps)	1
7-013-59	Ø6.0mm Trocar Tip Headed Reamer	1
7-013-60	Ø7.0mm Trocar Tip Headed Reamer	1
7-013-61	Ø8.0mm Trocar Tip Headed Reamer	1
7-013-62	Ø9.0mm Trocar Tip Headed Reamer	1
7-013-69	Awl/Tap for 4.75/5.5mm Gyrolock Screw-In Anchor	1
7-013-70	Tap for 4.75/5.5mm Gyrolock Screw-In Anchor	1
7-013-71	Tap for 3.5mm Gyrolock Screw-In Anchor	1
7-013-83	Ø2.5mm X 200mm Drill Bit	1
7-013-84	Ø3.2mm X 200mm Drill Bit	1
7-013-66	Instrument Tray A for Shoulder Arthroscopy Instrument Set	1
7-013-67	Instrument Tray B for Shoulder Arthroscopy Instrument Set	1
7-013-68	Instrument Tray C for Shoulder Arthroscopy Instrument Set	1
7-013-46	Container For Shoulder Arthroscopy Instrument Set	1
7-013-52	Lid For Container Shoulder Arthroscopy Instrument Set	1

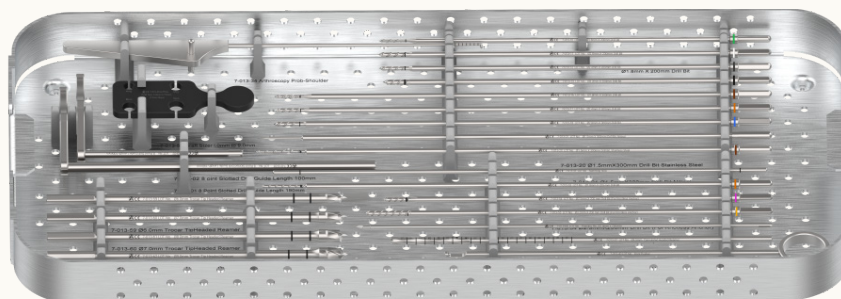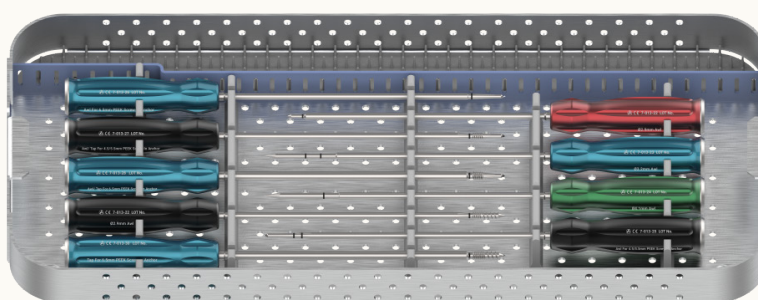7-013-66**Instrument Tray A for Shoulder Arthroscopy Instrument set****7-013-67****Instrument Tray B for Shoulder Arthroscopy Instrument set****7-013-68****Instrument Tray C for Shoulder Arthroscopy Instrument set**



7-013-46

Container For Shoulder Arthroscopy Instrument Set



7-013-52

Lid For Container Shoulder Arthroscopy Instrument Set

